


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Thalidomide

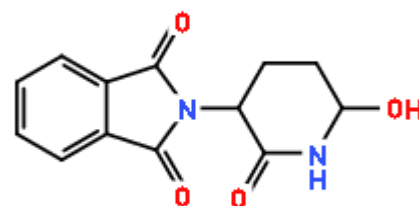

## CHROMATOGRAM

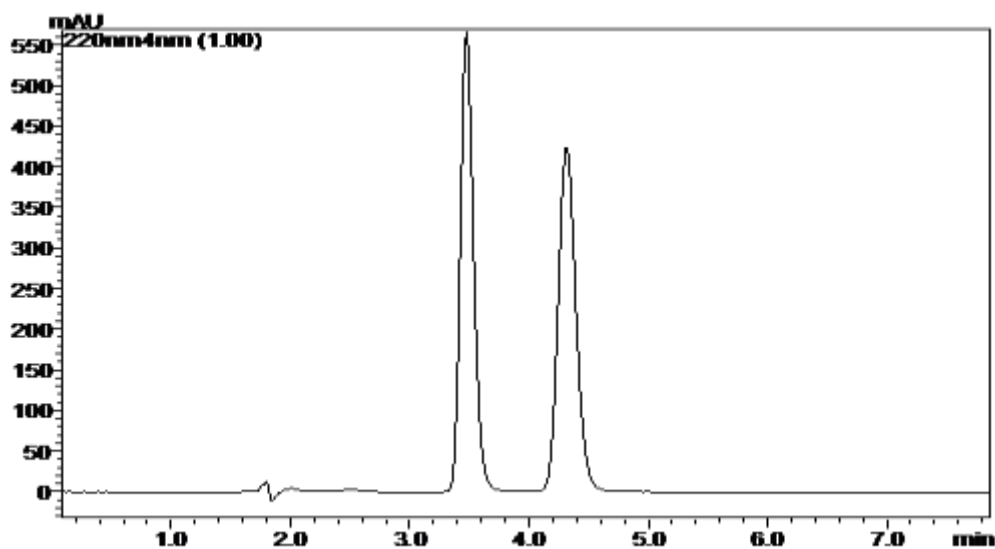

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® AD-3R		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	100mM ammonium acetate aq. / methanol / acetonitrile = 5 / 25 / 70		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	0.93	1 - 2	1.50
2	4.3	1.39		3.33