


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Propranolol

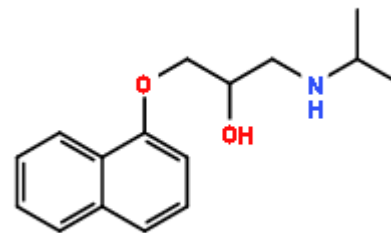

## CHROMATOGRAM

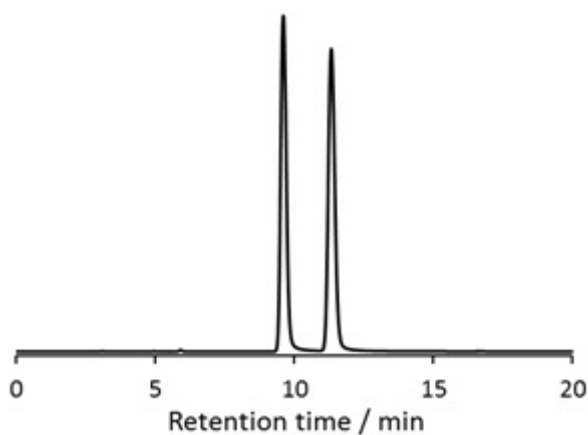

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	n-hexane / ethanol / diethylamine = 93 / 7 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.6	2.20	1 - 2	4.34
2	11.3	2.78		