


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Cetirizine

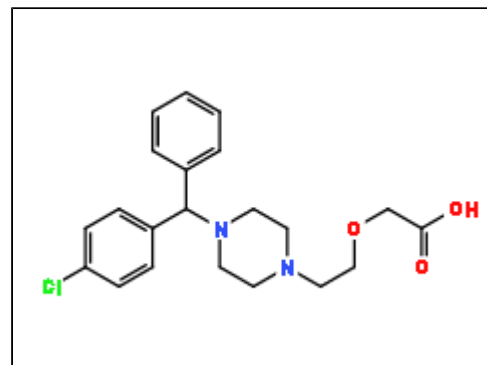

CHROMATOGRAM

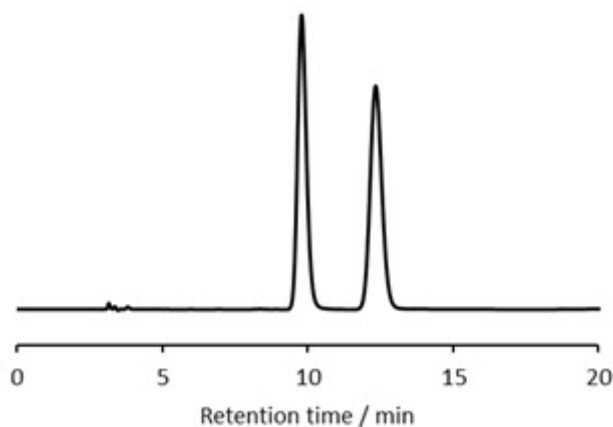

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	n-hexane / ethanol / acetic acid / diethylamine = 50 / 50 / 0.1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (μL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.8	2.26	1 - 2	3.87
2	12.3	3.11		