


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

4-Phenyl-2-oxazolidinone

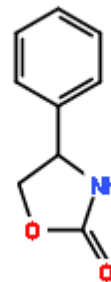

CHROMATOGRAM

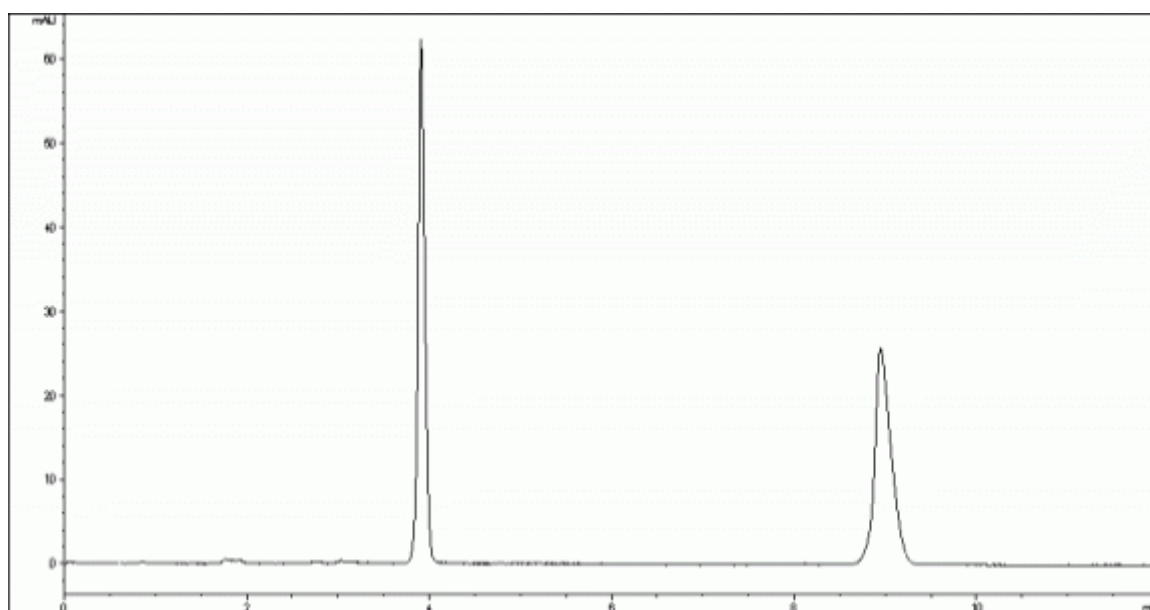

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	dichloromethane / ethanol = 98 / 2		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.02 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.9	1.30	1 - 2	22.94
2	9.0	4.26		