


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Frovatriptan

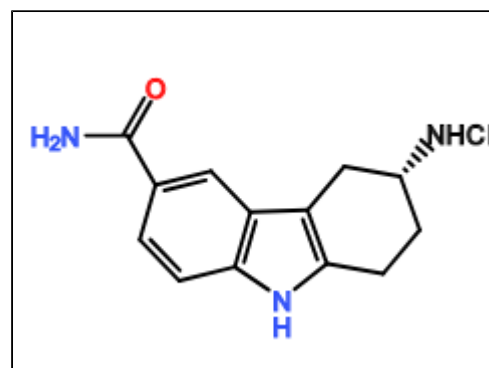

CHROMATOGRAM

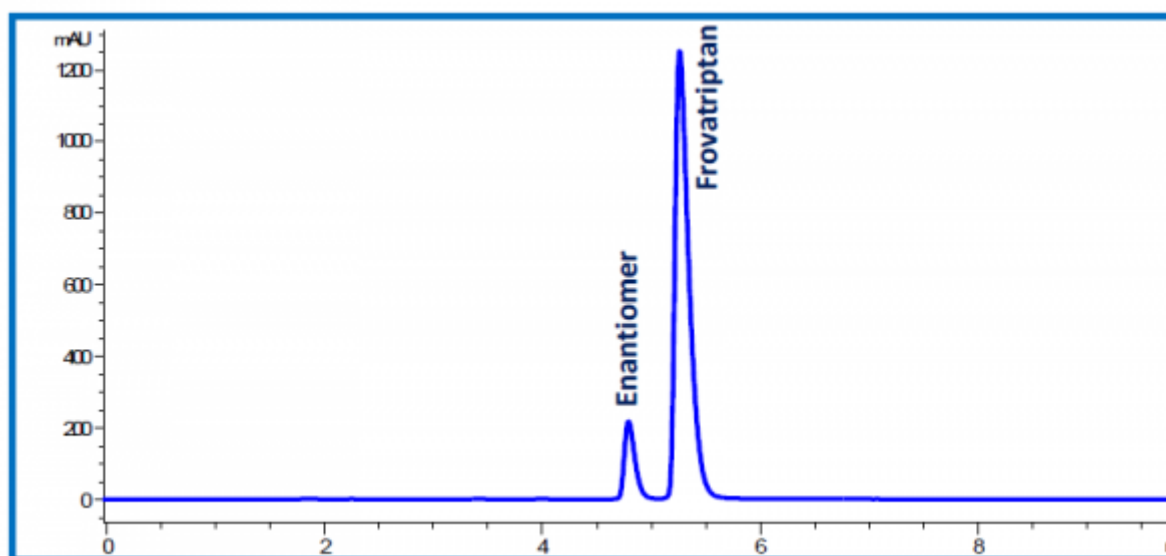

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	100mM ammonium bicarbonate aq. / acetonitrile = 10 / 90		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	0.010 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.8	1.82	1 - 2	2.30
2	5.3	2.12		