


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Aliskiren

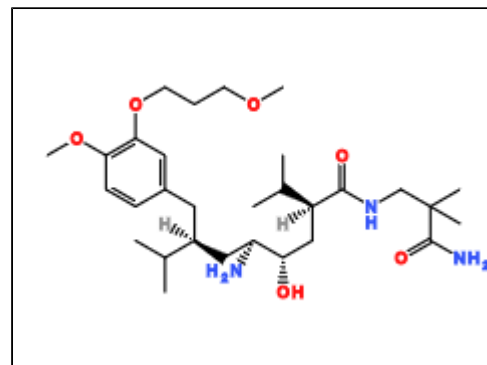

CHROMATOGRAM

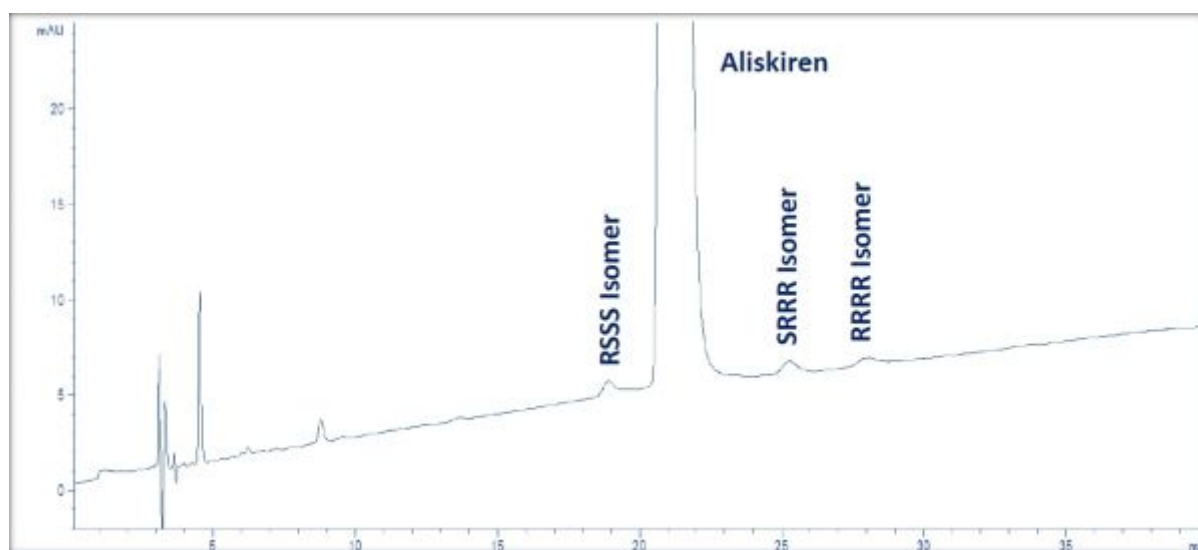

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IE-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	methyl tert-butyl ether / acetonitrile / 1,4-dioxane / ethylenediamine = 80 / 20 / 1 / 0.05		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	18.9		1 - 2	2.30
2	21.0		2 - 3	4.10
3	25.3		3 - 4	2.00
4	28.0			