


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Baclofen

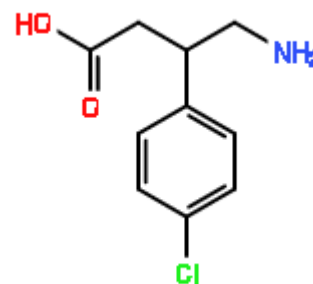

## CHROMATOGRAM

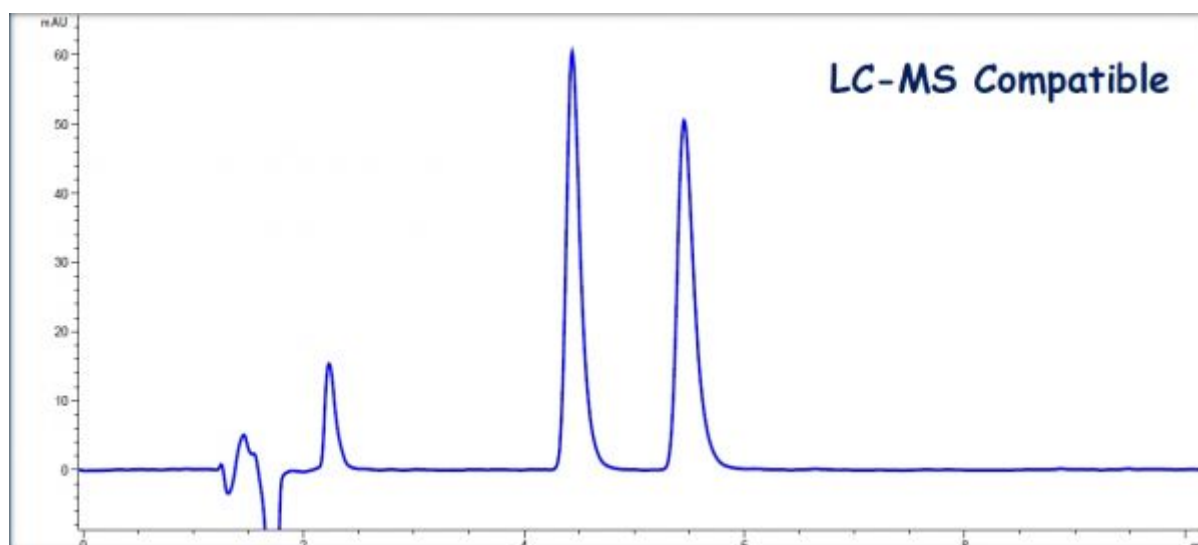

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 3 x 150 mm (3 µm)		
Mobile Phase Composition	methanol / acetonitrile / water / formic acid / ammonium formate (weight) = 750 / 200 / 50 / 0.2 / 0.32		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	0.007 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.4		1 - 2	4.02
2	5.5			