


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

barnidipine

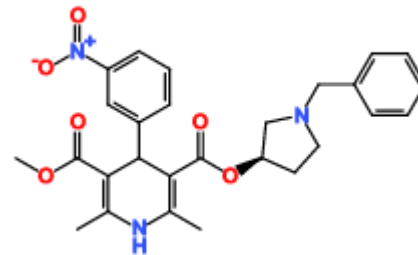

## CHROMATOGRAM

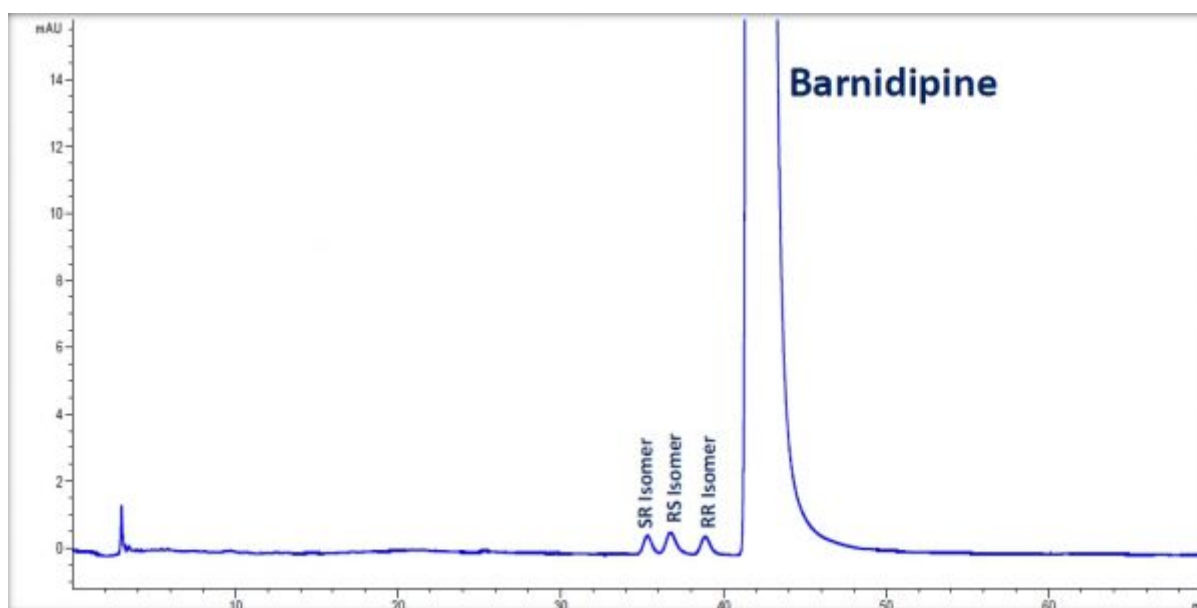

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	50mM ammonium bicarbonate aq. / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	35.3		1 - 2	1.37
2	36.7		2 - 3	1.97
3	38.9		3 - 4	2.23
4	41.8			