


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Bortezomib

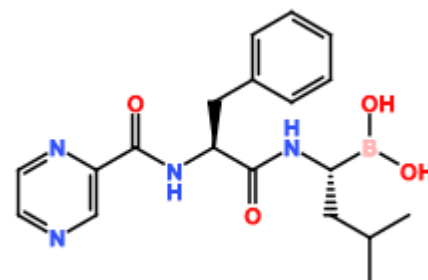

CHROMATOGRAM

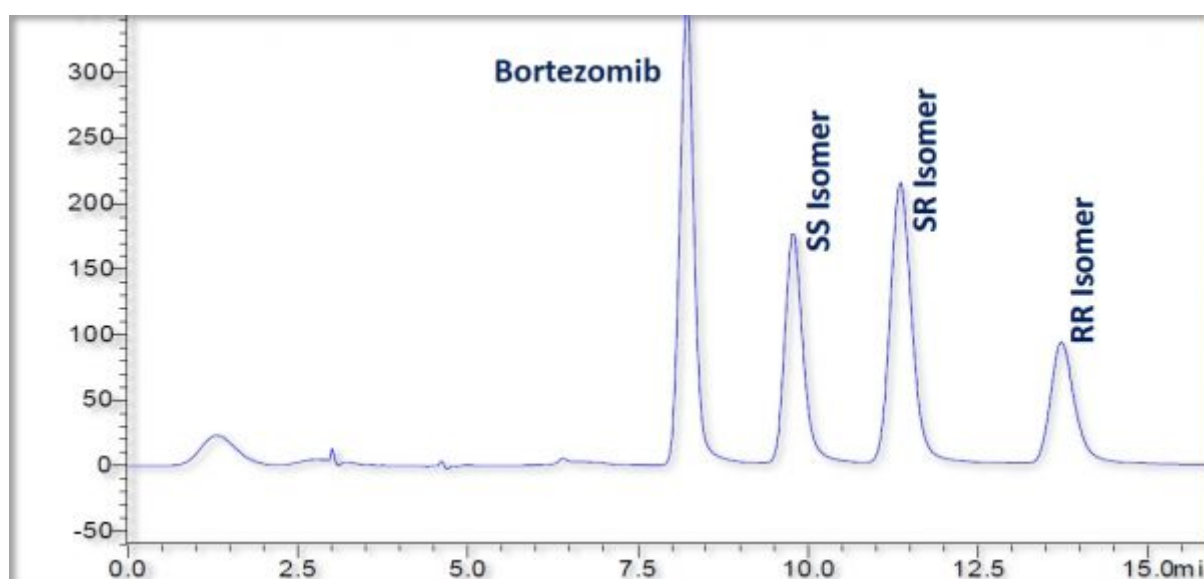

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® AD-3 4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-heptane / ethanol / 2-propanol / formic acid = 90 / 5 / 5 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 270 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.2		1 - 2	2.70
2	9.8		2 - 3	3.20
3	11.3		3 - 4	3.60
4	13.7			