


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Brivaracetam

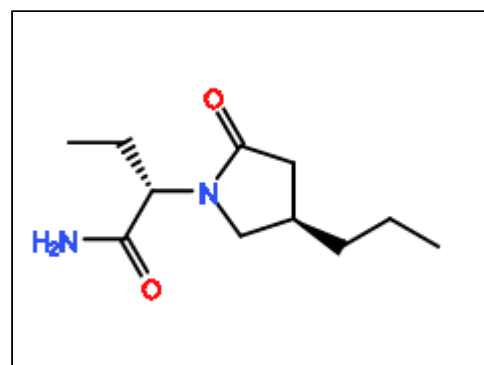

## CHROMATOGRAM

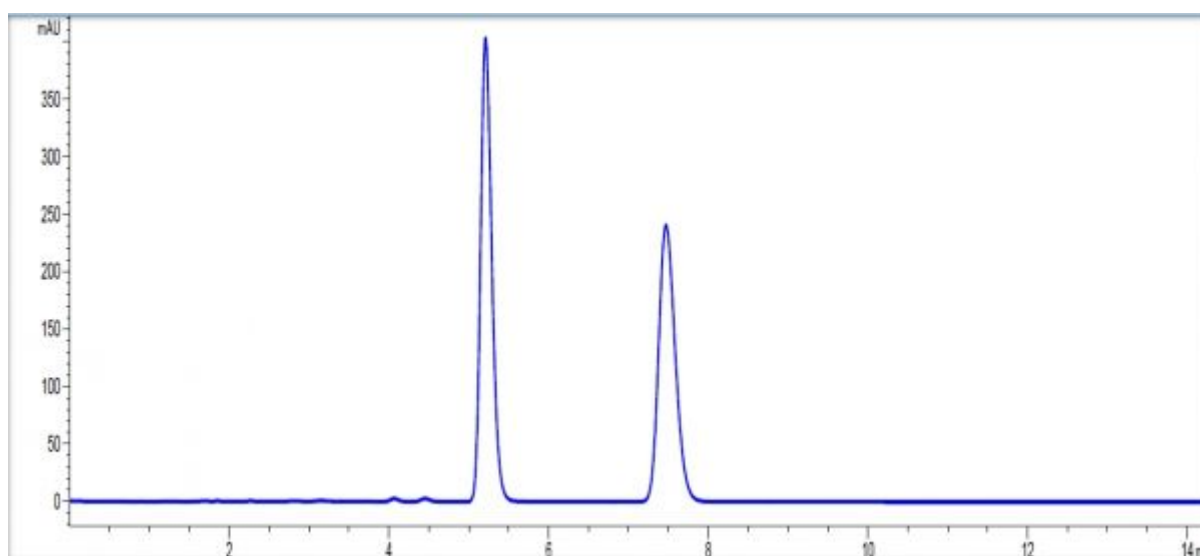

## ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® IG 4.6 x 150 mm (5 µm)		
Mobile Phase Composition	10mM ammonium bicarbonate aq. / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 215 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.2		1 - 2	7.10
2	7.4			