


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Eliglustat

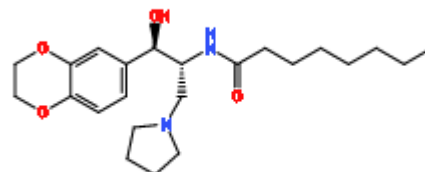

## CHROMATOGRAM

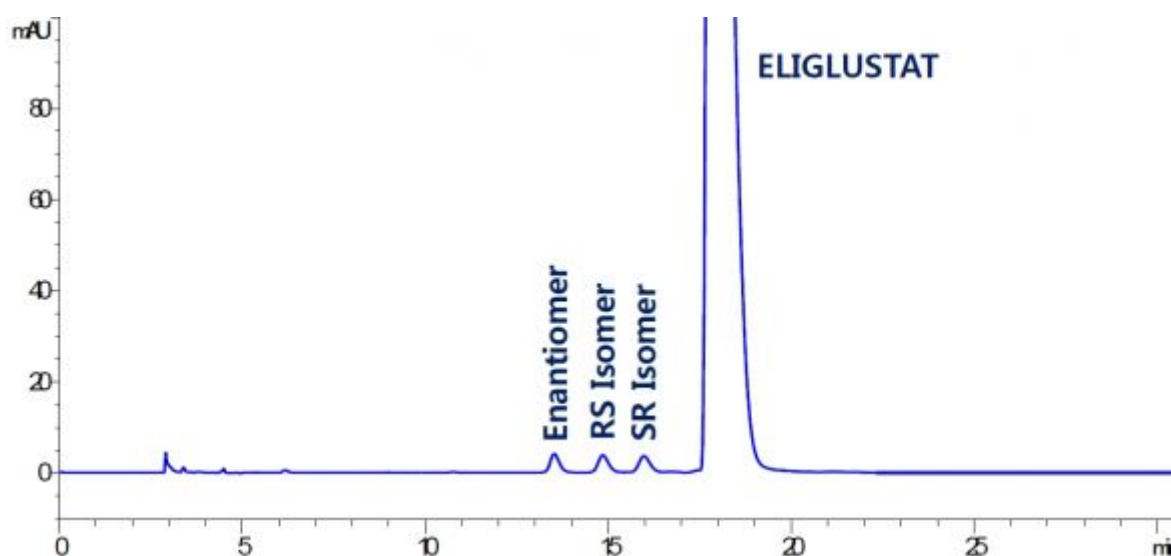

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IF-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / methanol / diethylamine = 90 / 5 / 5 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 280 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	13.5		1 - 2	2.50
2	14.9		2 - 3	2.10
3	16.0		3 - 4	2.30
4	17.8			