


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

5,5-Diphenyl-4-methyl-2-oxazolidinone

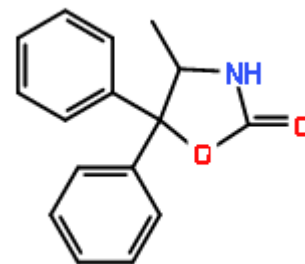

CHROMATOGRAM

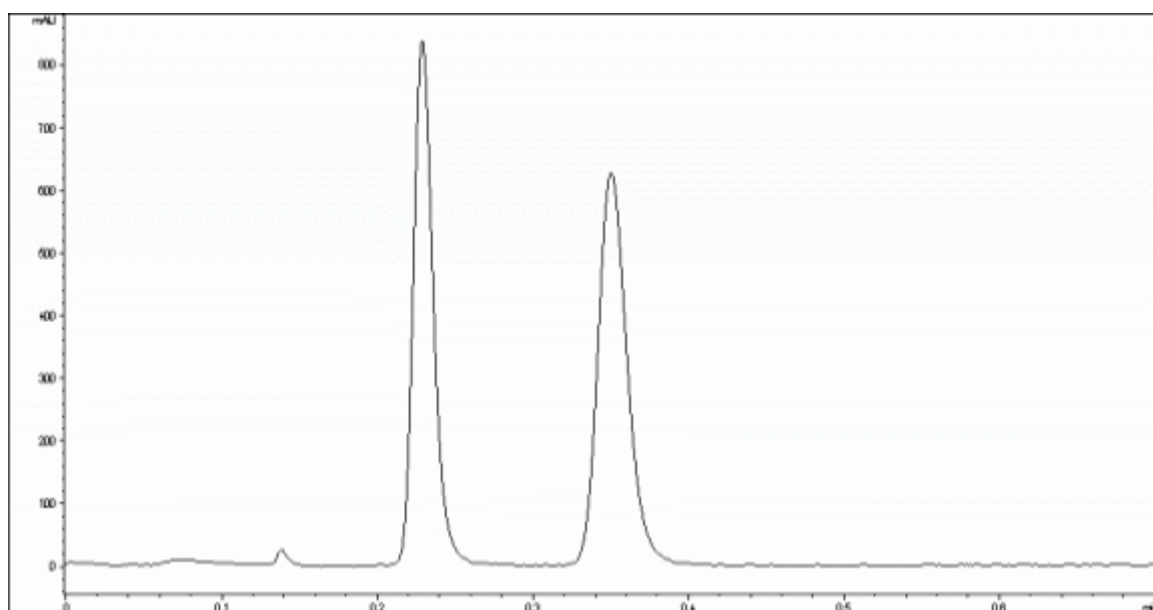

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	methyl tert-butyl ether / ethanol = 95 / 5		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.84	1 - 2	4.28
2	0.4	1.82		