


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Cinacalcet

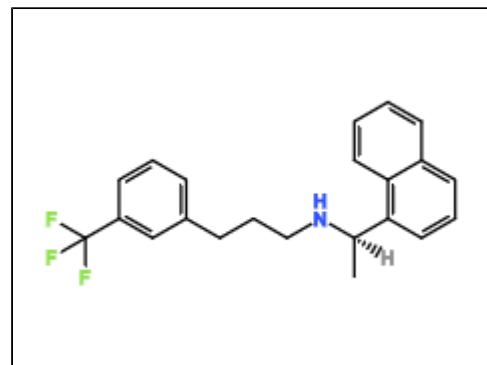

CHROMATOGRAM

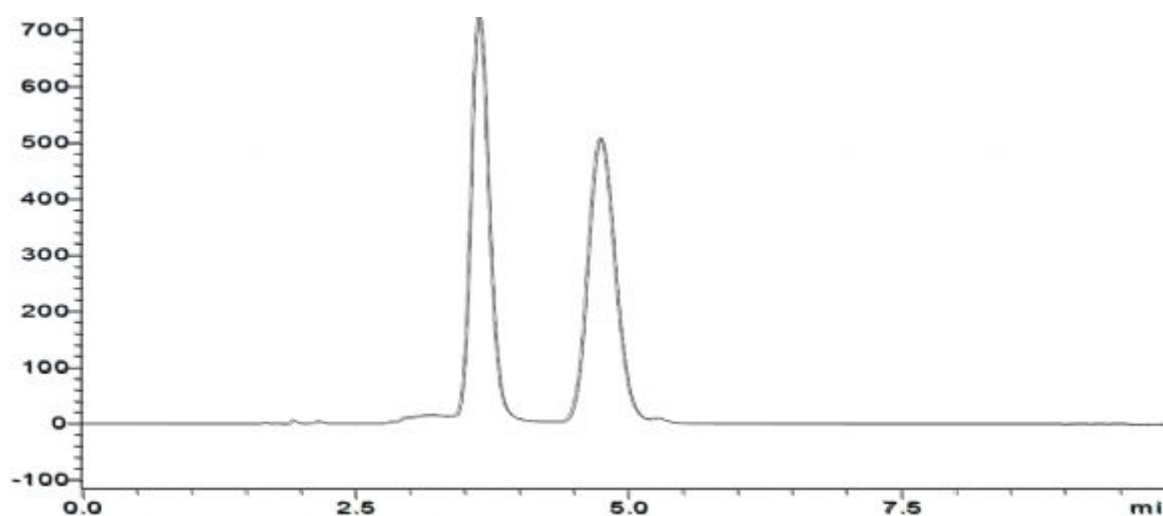

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® AD-3R 4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / methanol / ammonium bicarbonate (weight) = 100 / 900 / 0.4		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.01 (mg)	Sample Concentration	
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.6		1 - 2	2.90
2	4.7			