


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Dexrazoxane

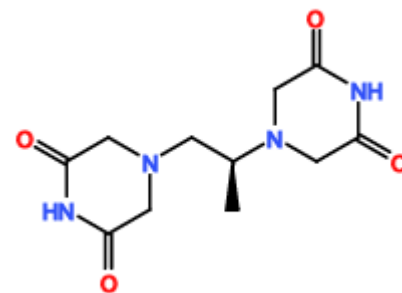

## CHROMATOGRAM

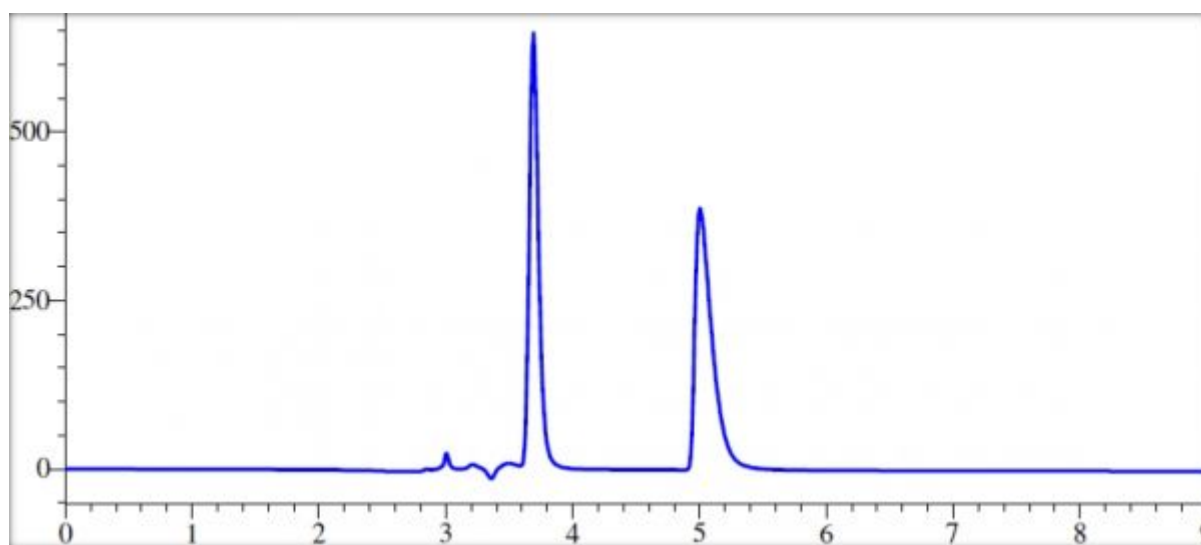

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ID		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	methanol / acetonitrile / diethylamine = 60 / 40 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 235 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.7		1 - 2	6.47
2	5.0			