


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Efavirenz

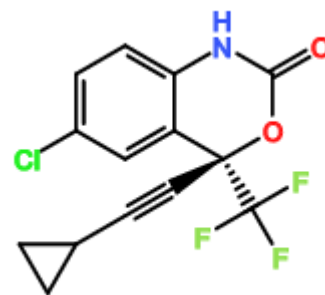

CHROMATOGRAM

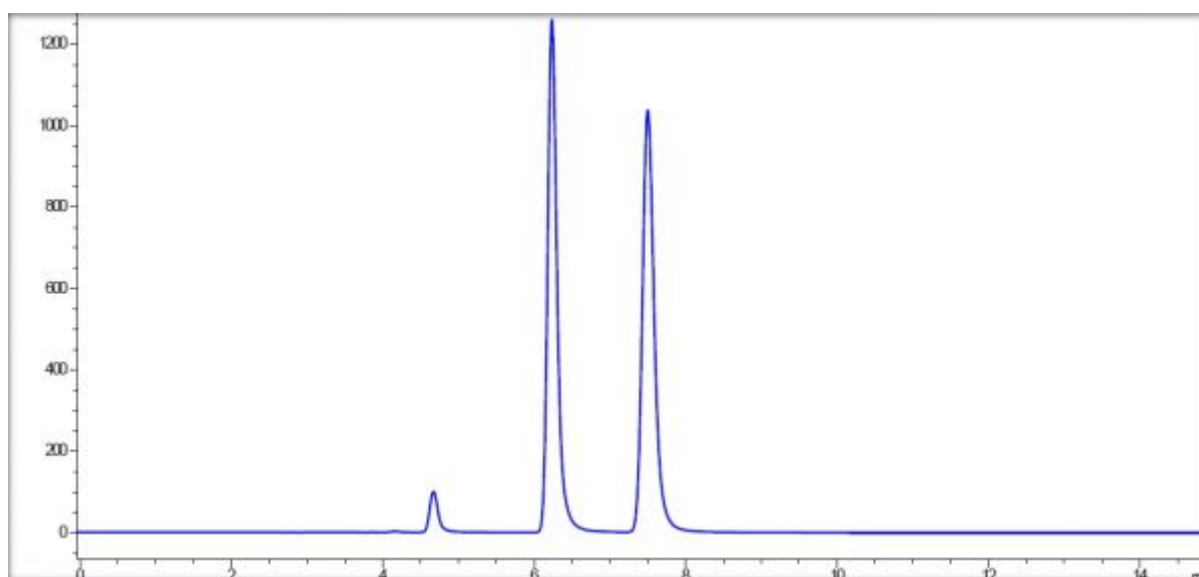

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	water / methanol / ammonium bicarbonate (weight) = 200 / 800 / 0.4		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 246 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.2		1 - 2	5.20
2	7.5			