


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ivabradine

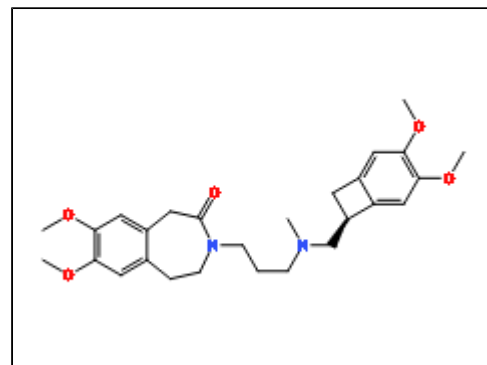

## CHROMATOGRAM

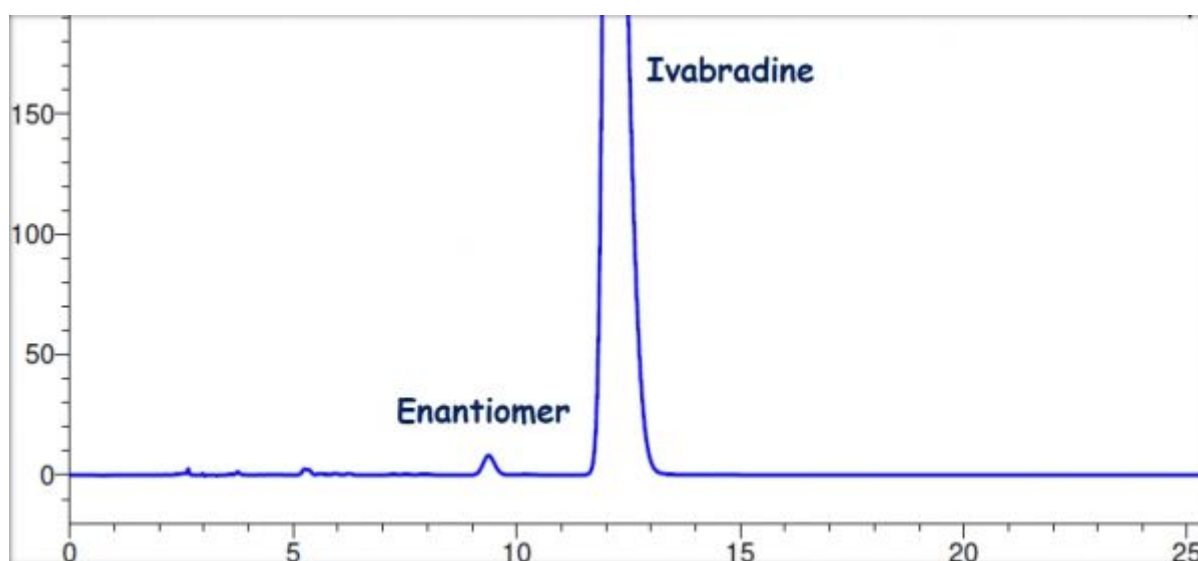

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	methanol / acetonitrile / diethylamine = 50 / 50 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (μL)		Sample Concentration	
Detection	UV-VIS 285 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.4		1 - 2	4.00
2	12.1			