


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Labetalol

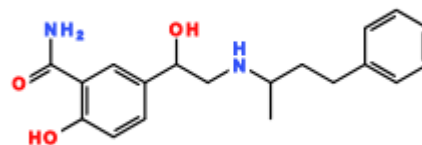

CHROMATOGRAM

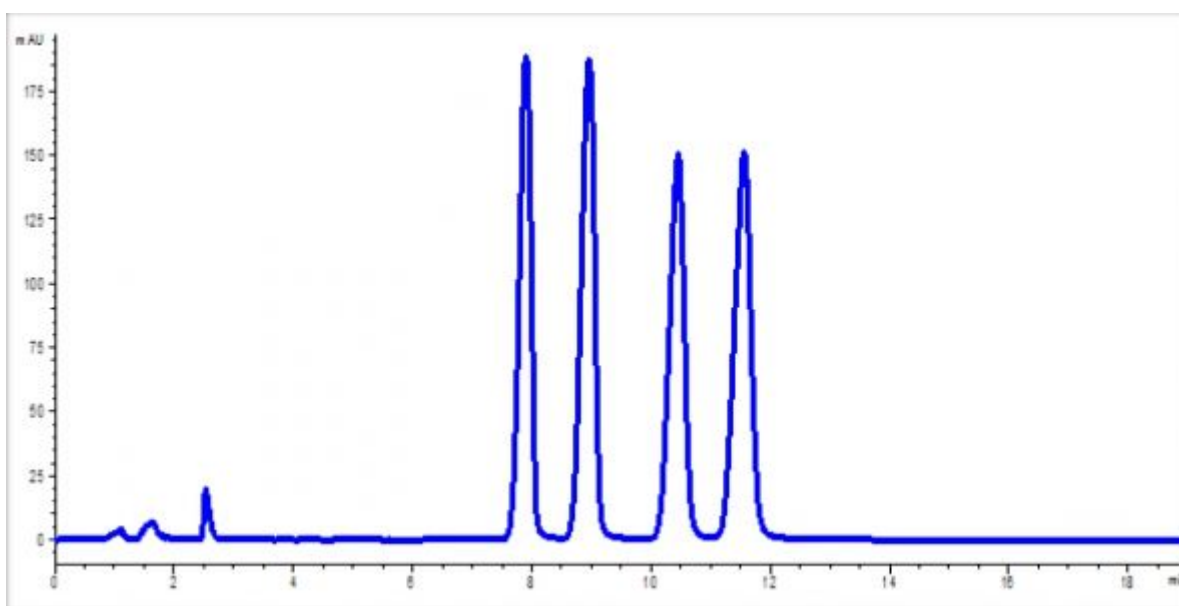

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(+)		
Column Size (Particle Size)	4 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / tetrahydrofuran / water / triethylamine = 490 / 490 / 20 / 1.74		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 310 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.9		1 - 2	2.70
2	9.0		2 - 3	3.40
3	10.5		3 - 4	2.20
4	11.6			