


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Levothyroxine

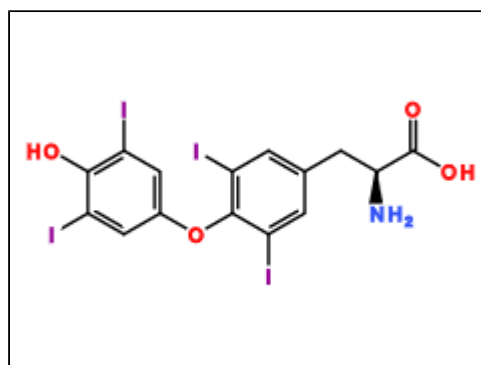

## CHROMATOGRAM

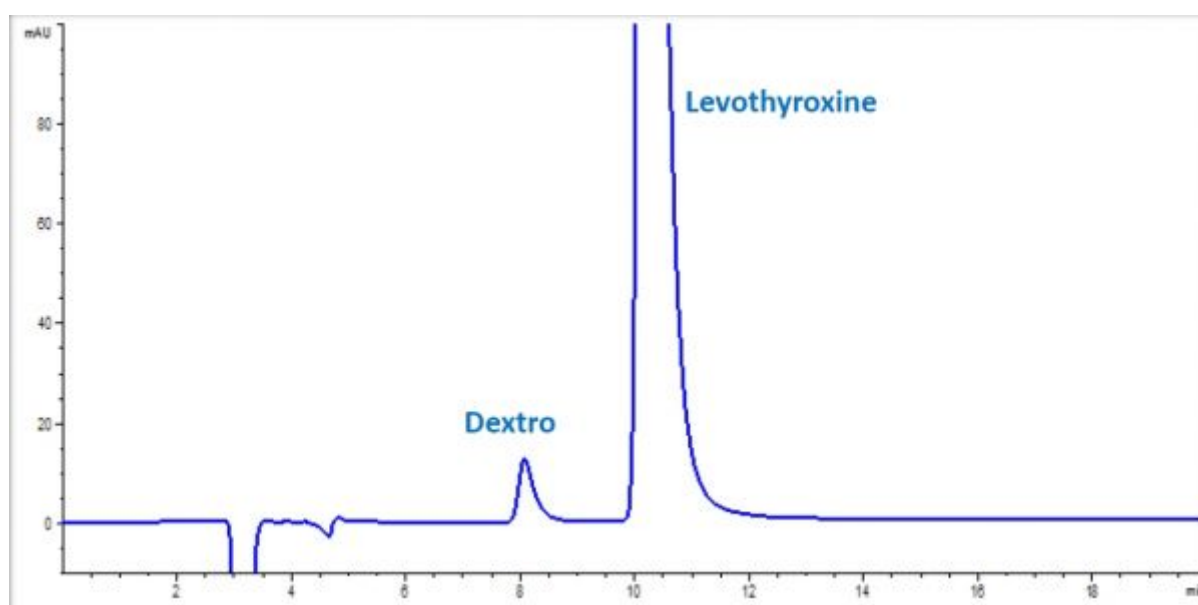

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® ZWIX(-)		
Column Size (Particle Size)	3 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / tetrahydrofuran / water / formic acid / diethylamine = 450 / 500 / 50 / 1.9 / 2.6		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.0		1 - 2	3.91
2	10.2			