


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Linezolid

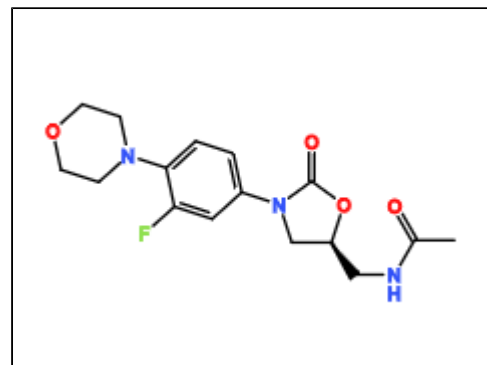

CHROMATOGRAM

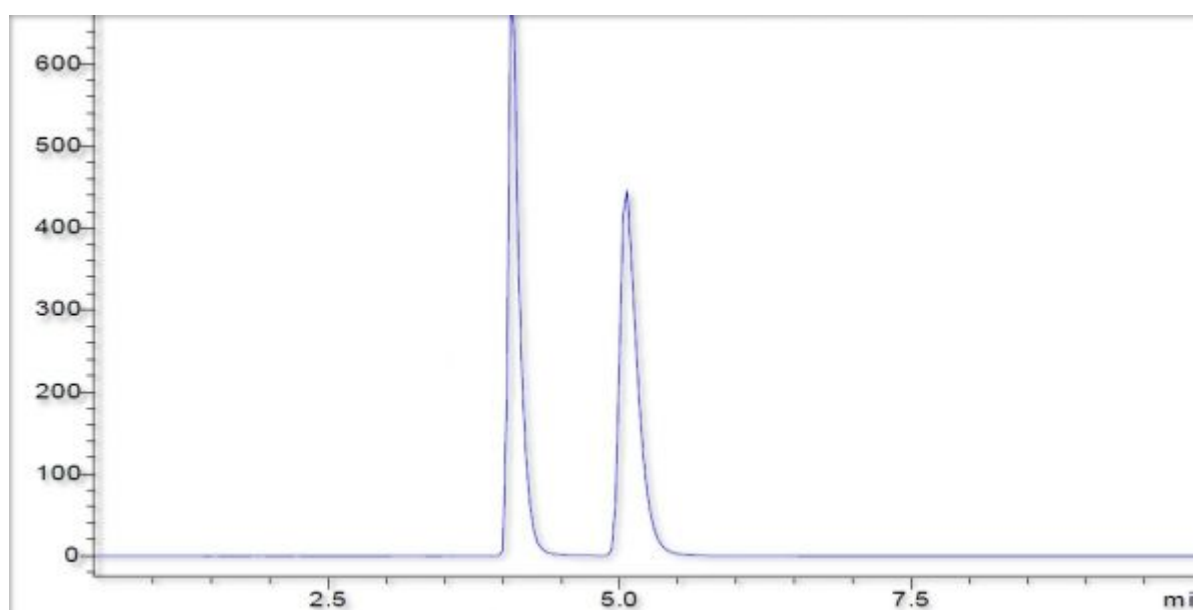

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	acetonitrile / methanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 260 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.1		1 - 2	3.80
2	5.0			