


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

alpha-Methyl-DL-tyrosine

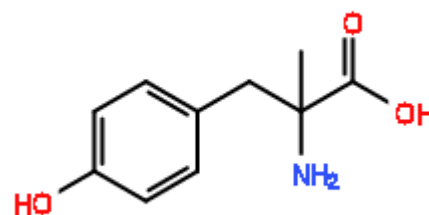

CHROMATOGRAM

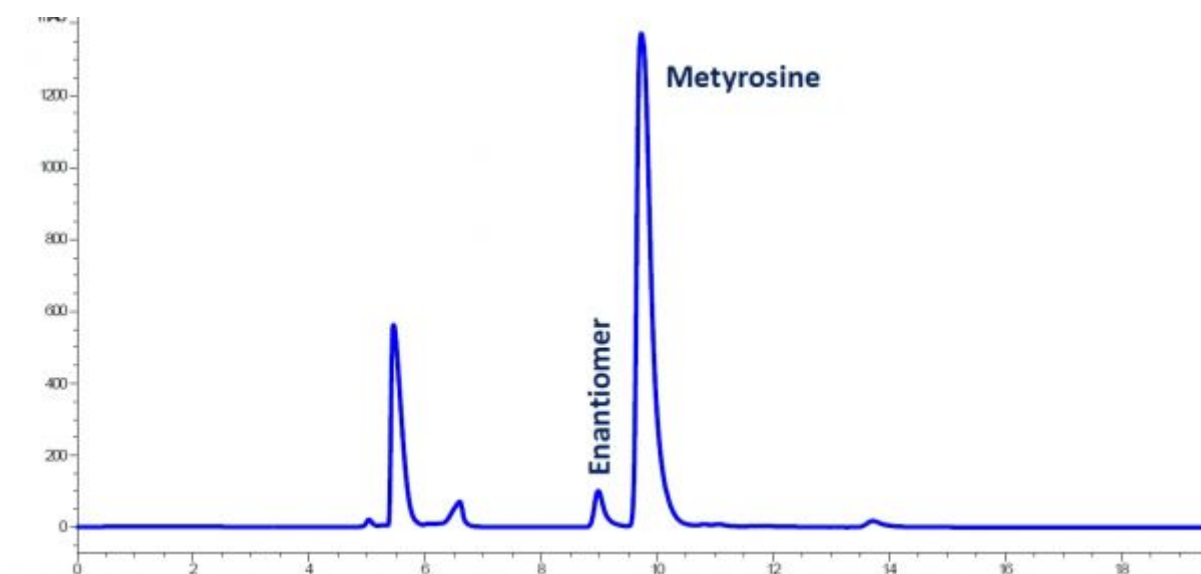

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CHIRALPAK® ZWIX(+) 4 x 250 mm (3 µm)		
Mobile Phase Composition	methanol / water / formic acid / diethylamine = 980 / 20 / 1.9 / 2.6		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 224 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.0		1 - 2	2.10
2	9.7			