


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ranolazine

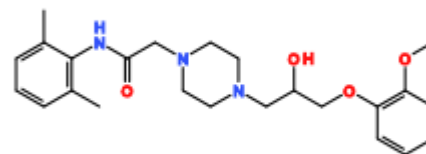

## CHROMATOGRAM

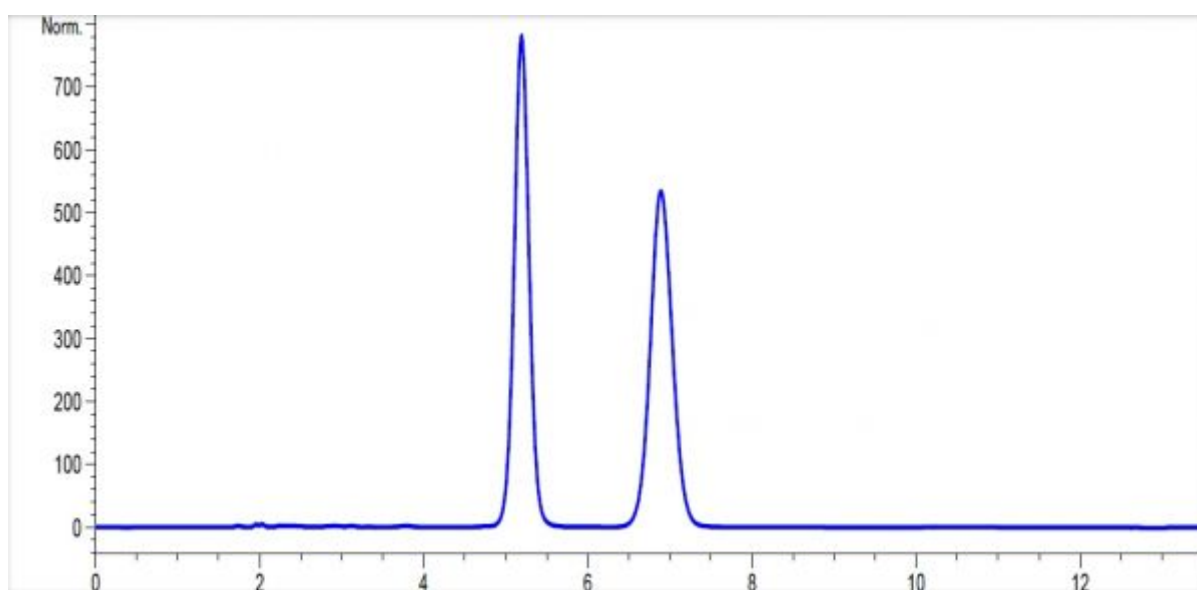

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / methanol / ammonium bicarbonate (weight) = 50 / 950 / 0.79		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.2		1 - 2	4.12
2	6.9			