


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Saxagliptin

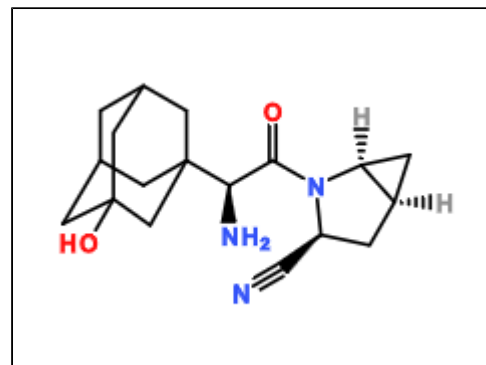

CHROMATOGRAM

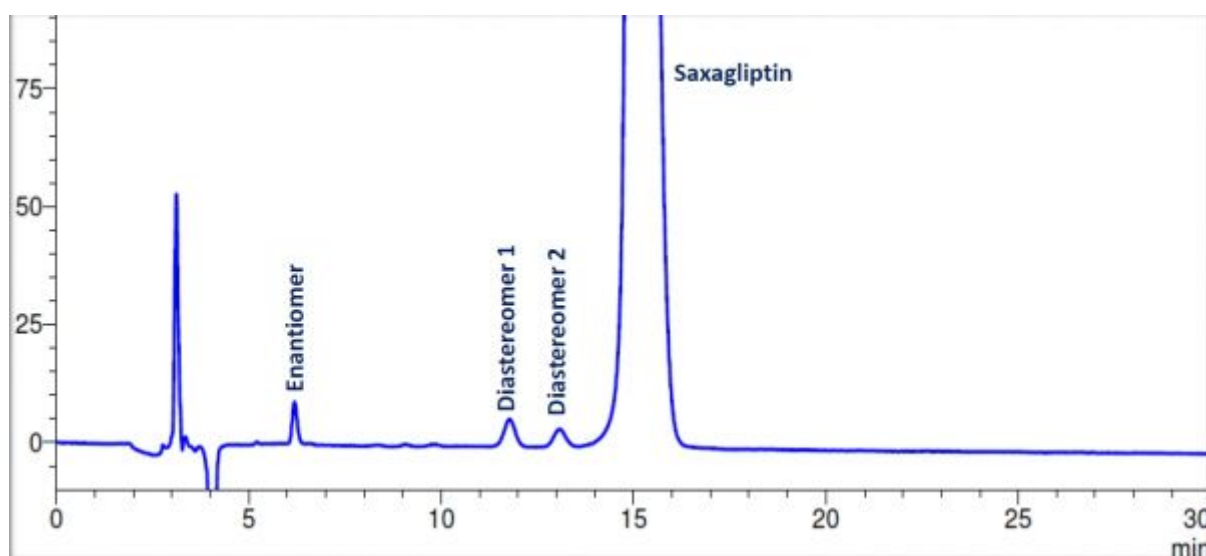

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	ethanol / acetonitrile / isopropylamine = 20 / 80 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.2		1 - 2	14.10
2	11.8		2 - 3	2.30
3	13.1		3 - 4	2.50
4	15.1			