


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Solifenacin (containing two pairs of enantiomers)

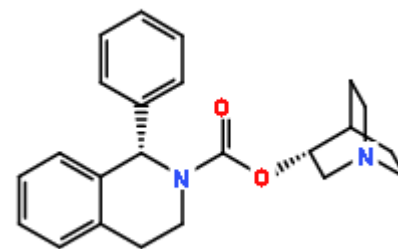

CHROMATOGRAM

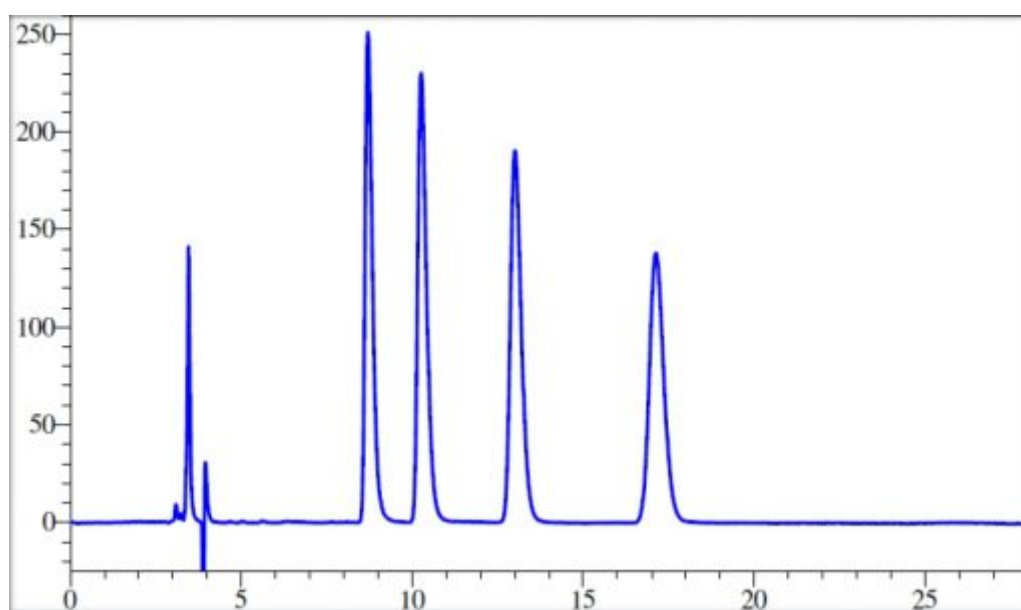

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / 2-propanol / diethylamine = 60 / 15 / 25 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	30
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.7		1 - 2	3.50
2	10.3		2 - 3	5.10
3	13.0		3 - 4	6.00
4	17.1			