


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Phenylethylamine

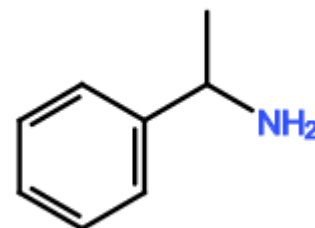

CHROMATOGRAM

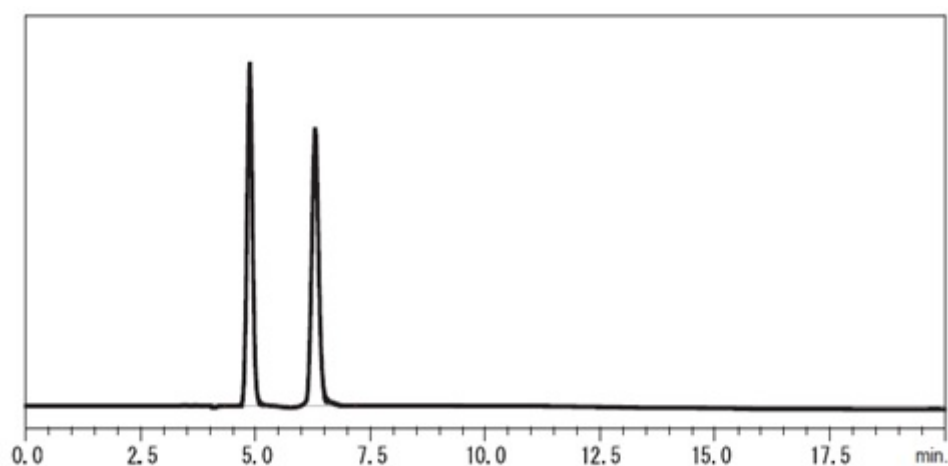

ANALYTICAL REPORT

Column Trade Name Column Size (Particle Size)	CROWNPAK® CR-I(+) 3 x 150 mm (5 µm)		
Mobile Phase Composition	perchloric acid aq. pH1.5 / acetonitrile = 50 / 50		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.2	Temperature (°C)	25
Injection Volume (µL)	0.001 (mg)	Sample Concentration	
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1			1 - 2	