


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

9,10-Dihydro-2-methyl-4H-benzo-[5,6] cyclohept[1,2-d]oxazol-4-ol

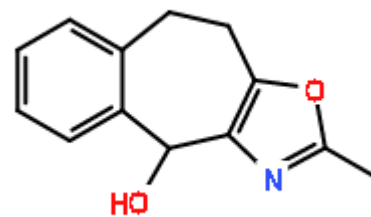

CHROMATOGRAM

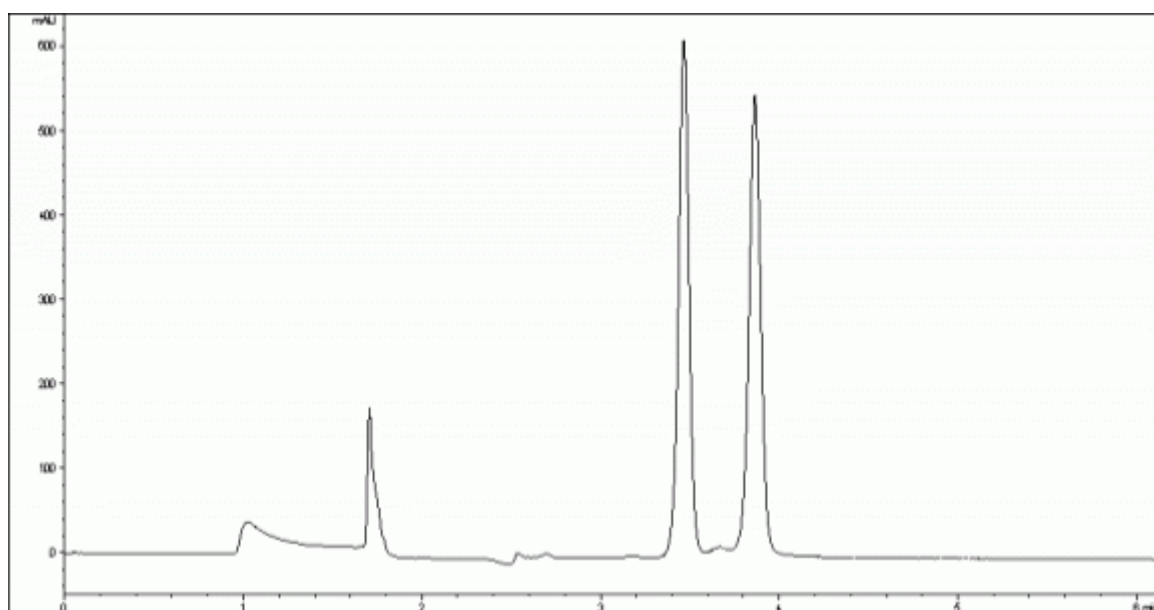

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / ethanol / ethanalamine = 50 / 50 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.5	1.04	1 - 2	3.59
2	3.9	1.27		