


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

alpha-(2,4-Dichlorophenyl)-1H-imidazole-1-ethanol

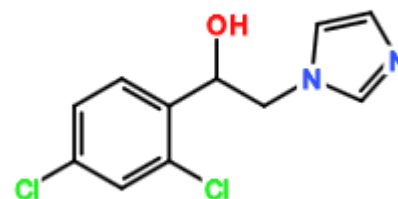

CHROMATOGRAM

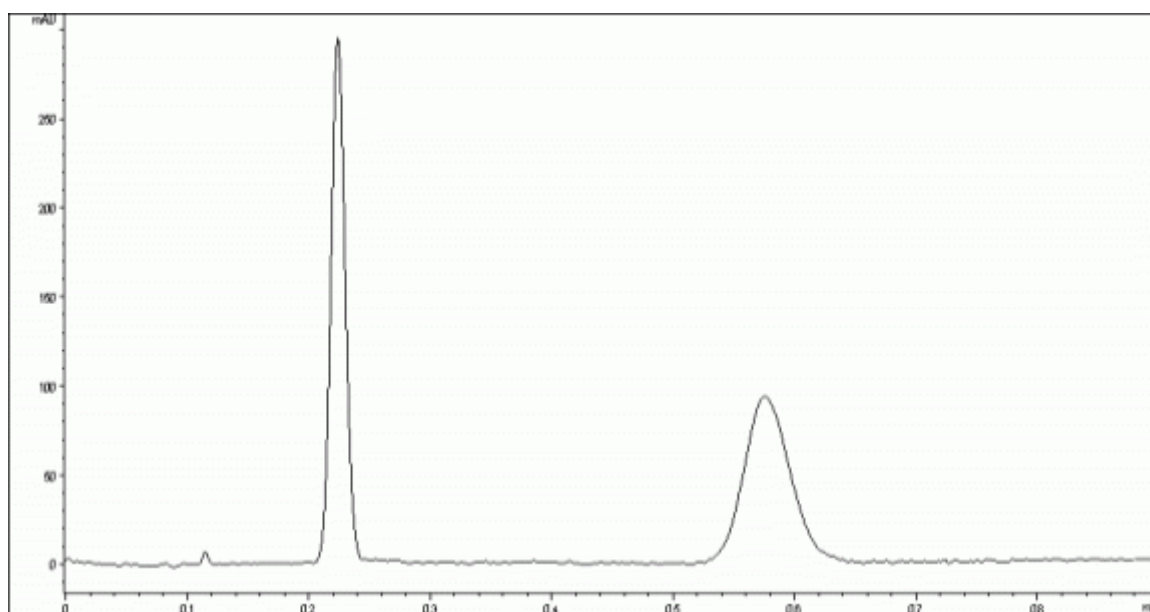

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 50 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / ethanolamine = 80 / 20 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	5	Temperature (°C)	25
Injection Volume (µL)	0.002 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	0.2	0.95	1 - 2	7.66
2	0.6	4.01		