


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Hydroxychloroquine Sulfate

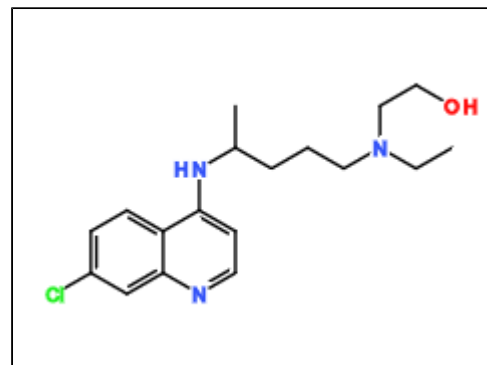

## CHROMATOGRAM

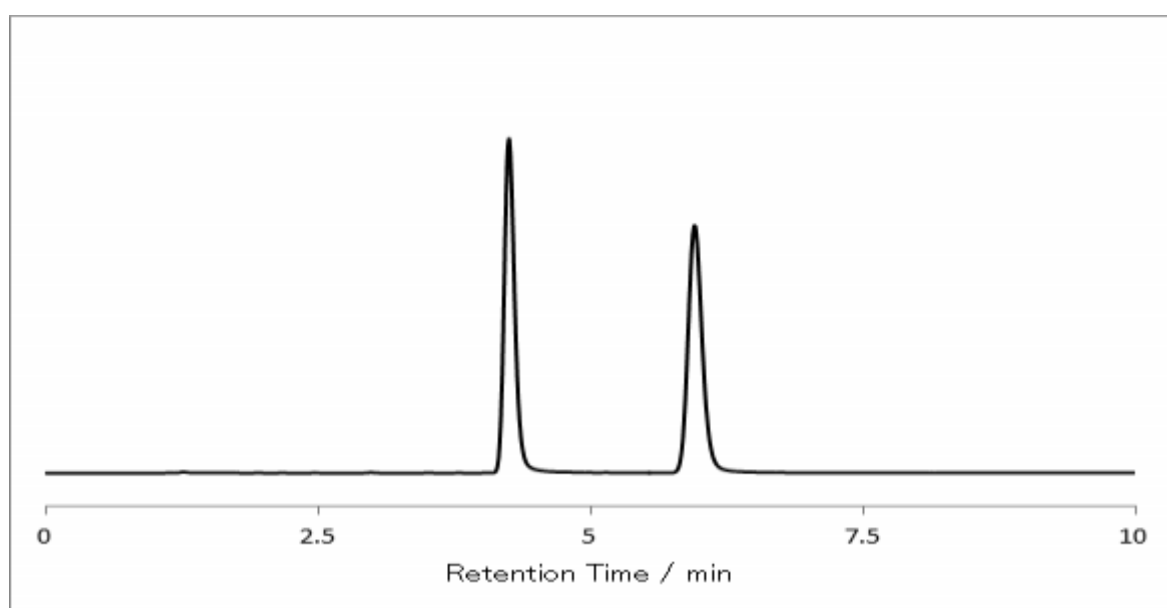

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	acetonitrile / methanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Polar Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	40
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 328 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.2	1.50	1 - 2	7.87
2	6.0	2.50		