


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Atorvastatin (racemic)

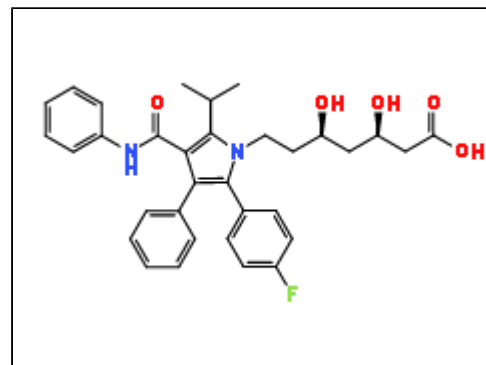

CHROMATOGRAM

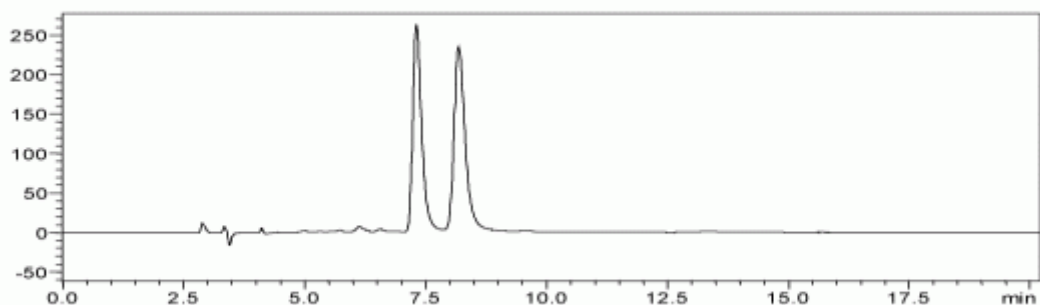

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / trifluoroacetic acid = 85 / 15 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 245 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.3	1.44	1 - 2	2.34
2	8.2	1.73		