


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Warfarin

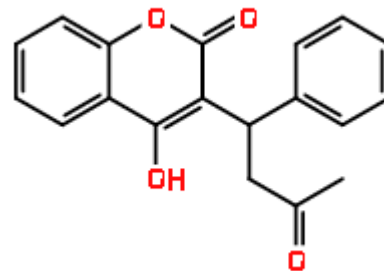

## CHROMATOGRAM

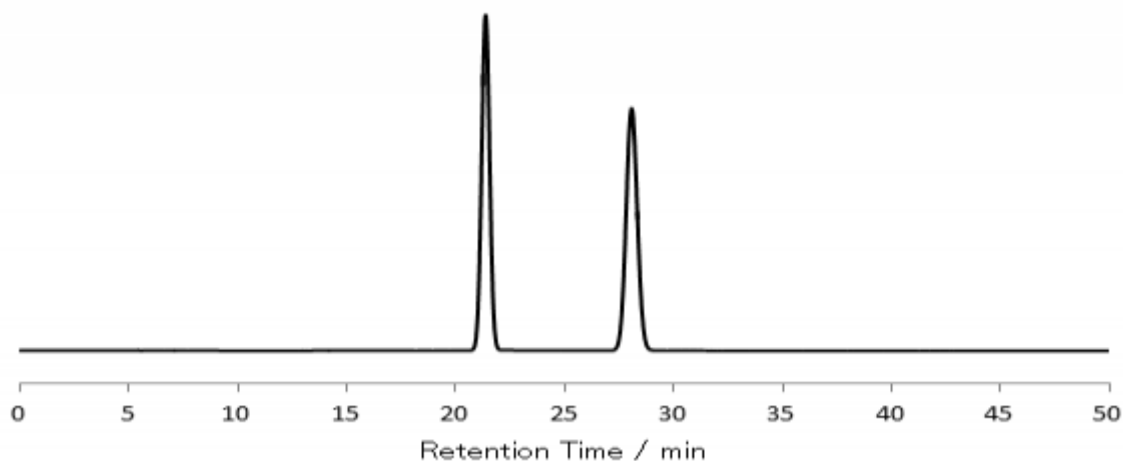

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB N-5		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	H3PO4 aq. pH2.0 / acetonitrile = 90 / 10		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.5	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	21.4	2.60	1 - 2	8.03
2	28.1	3.70		