

## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Phenyl vinyl sulfoxide



## CHROMATOGRAM

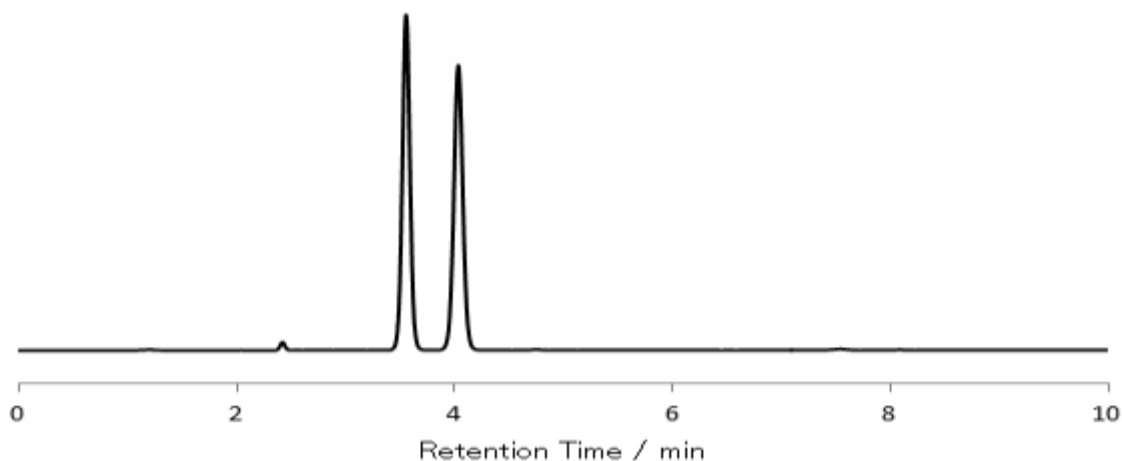

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	water / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.003 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.6	1.11	1 - 2	3.41
2	4.0	1.40		