


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

1-(1-Naphthyl)-ethanol

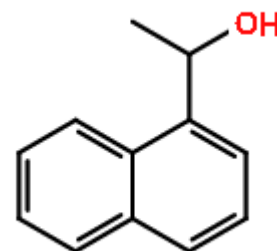

## CHROMATOGRAM

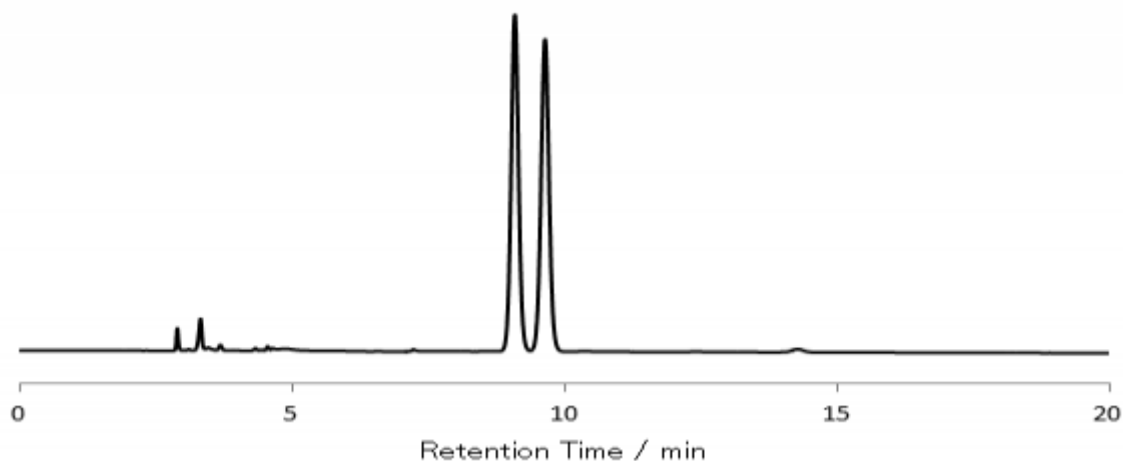

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / ethanol = 85 / 15 / 1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	9.1	2.20	1 - 2	2.14
2	9.6	2.30		