


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

gamma-Phenyl-gamma-butyrolactone

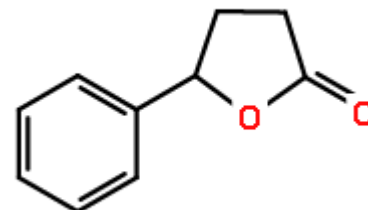

CHROMATOGRAM

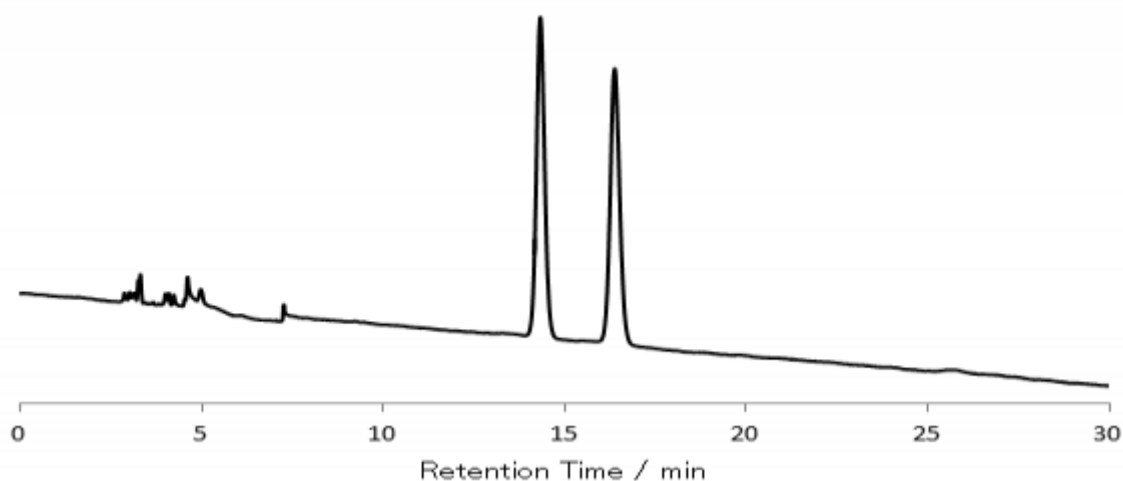

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / ethanol = 85 / 15 / 1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.3	4.00	1 - 2	4.66
2	16.4	4.70		