


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sulpiride

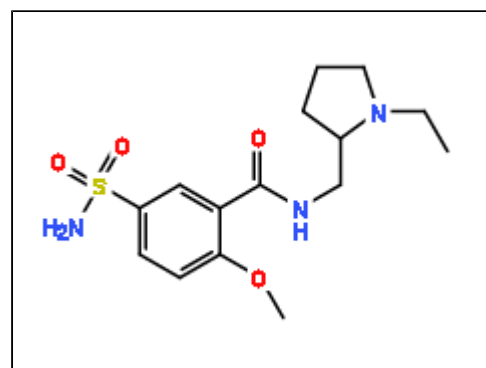

CHROMATOGRAM

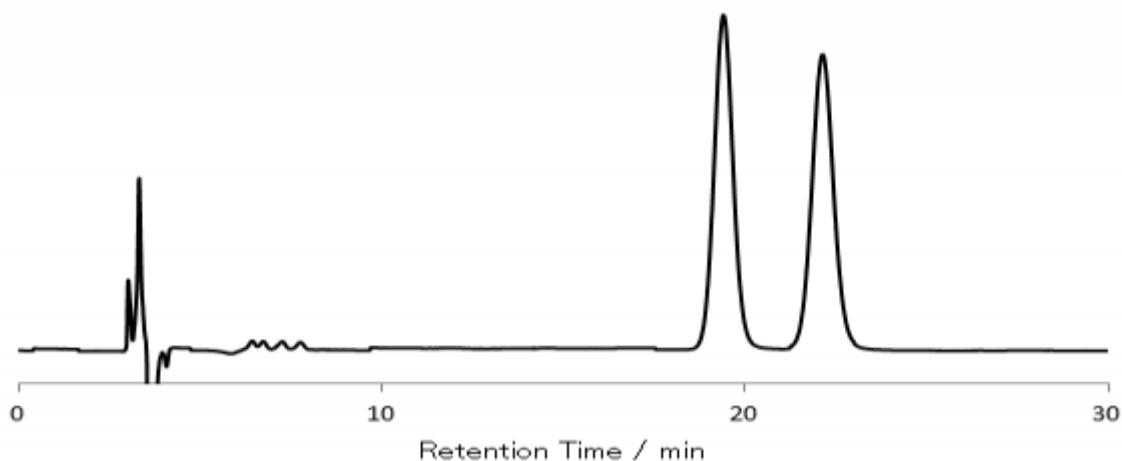

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / tetrahydrofuran / diethylamine = 60 / 40 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	19.4	5.70	1 - 2 1.17	2.67
2	22.1	6.60		