


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Ofloxacin Ethyl ester

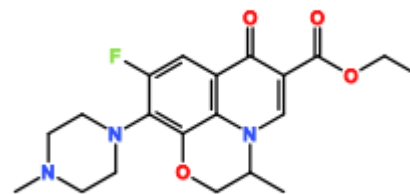

CHROMATOGRAM

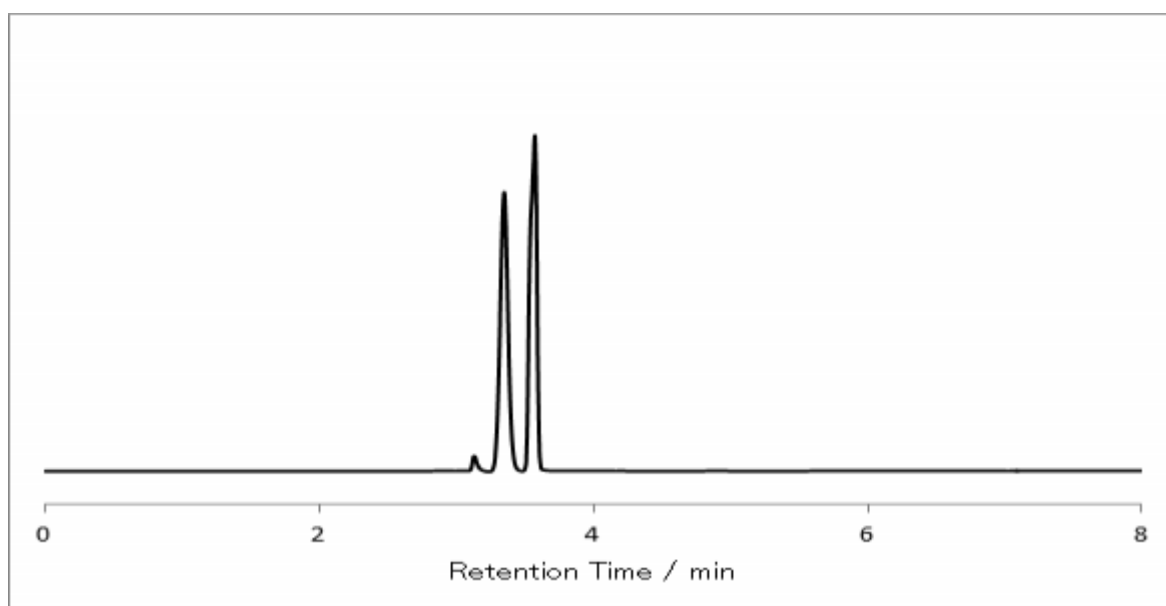

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / ethanol / diethylamine = 10 / 90 / 1 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 300 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.4	0.16	1 - 2	1.48
2	3.6	0.24		2.12