


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Sulpiride

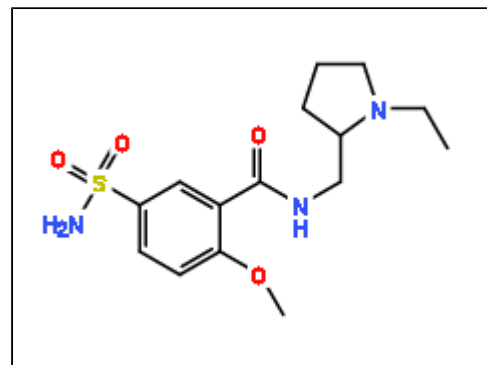

CHROMATOGRAM

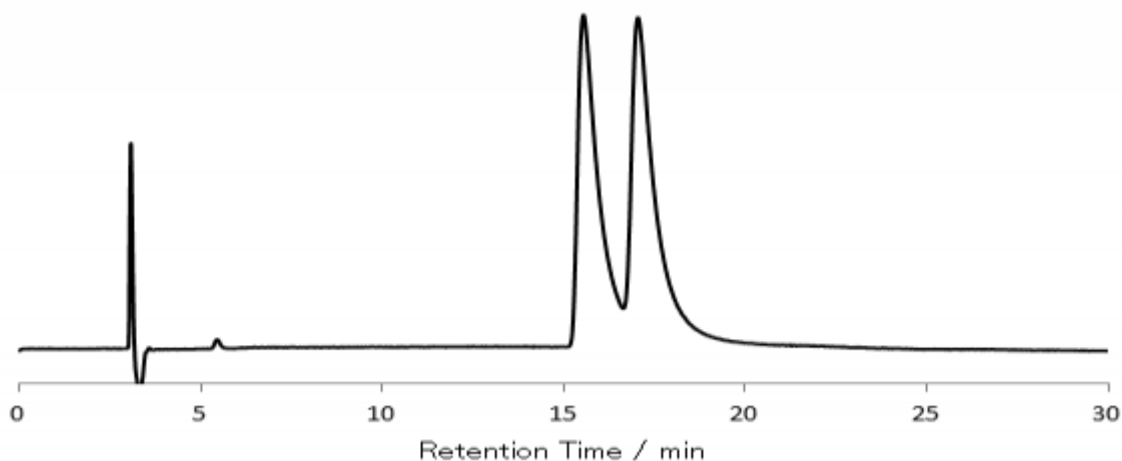

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethyl acetate / diethylamine = 70 / 30 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 284 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	15.5	4.18	1 - 2 1.12	1.50
2	17.0	4.68		