


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Meclizine

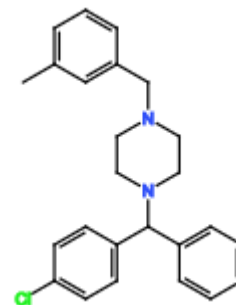

CHROMATOGRAM

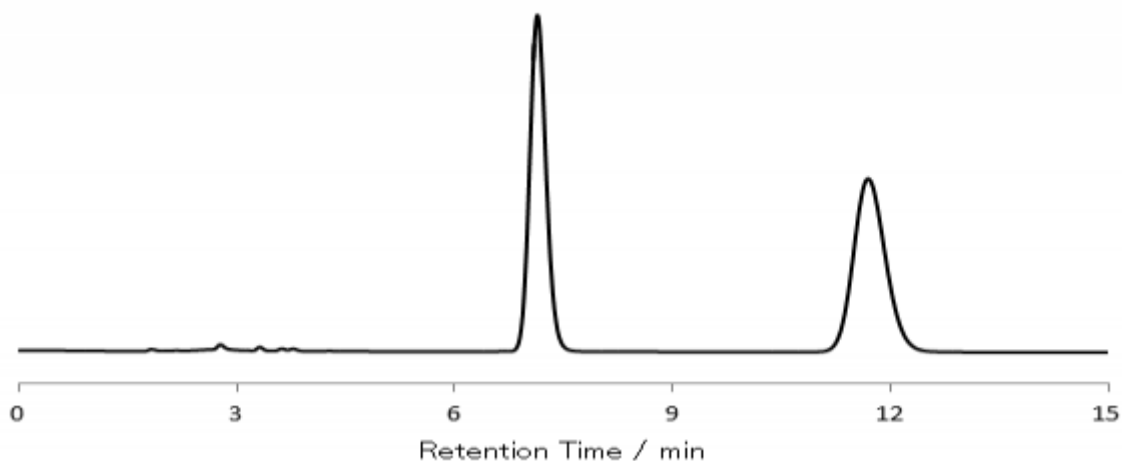

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 150 mm (5 µm)		
Mobile Phase Composition	20mM NH ₄ HCO ₃ aq. pH9.0 / acetonitrile = 20 / 80		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.1	2.57	1 - 2	7.42
2	11.7	4.84		