


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Ethotoin

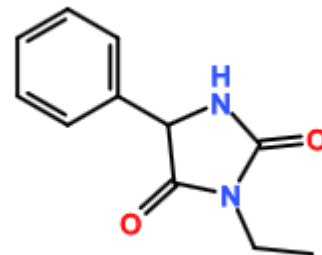

## CHROMATOGRAM

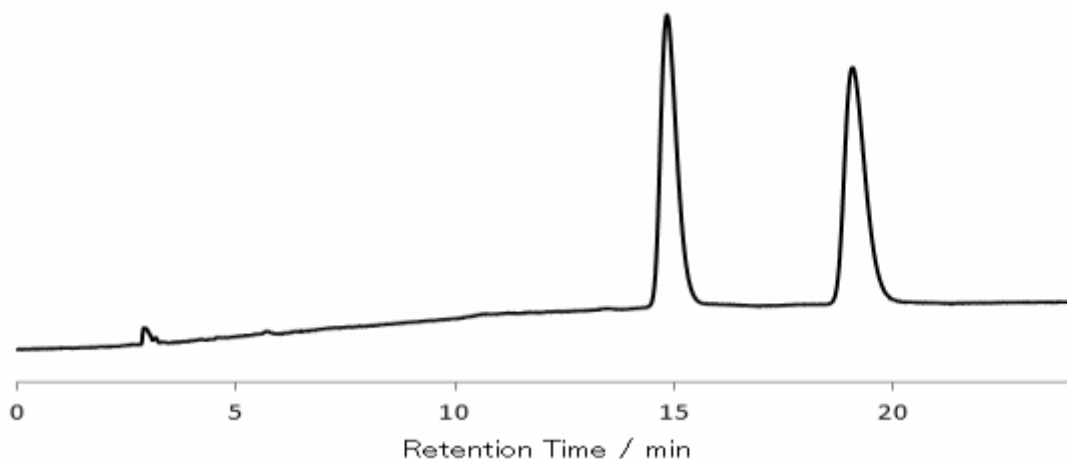

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol = 85 / 15		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 200 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	14.8	3.94	1 - 2	5.43
2	19.1	5.36		