


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Propranolol

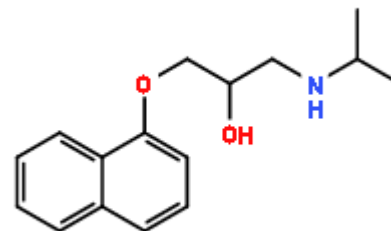

CHROMATOGRAM

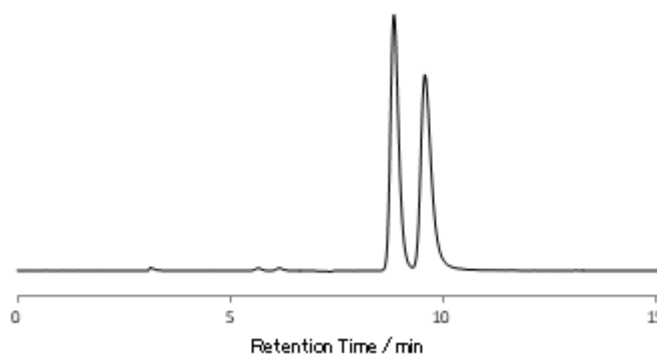

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IJ		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / ethanol / ethylenediamine = 95 / 5 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 290 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.9	1.95	1 - 2	1.83
2	9.6	2.19		