


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

4-Benzoyloxy-2-azetidinone

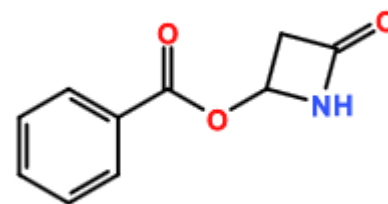

CHROMATOGRAM

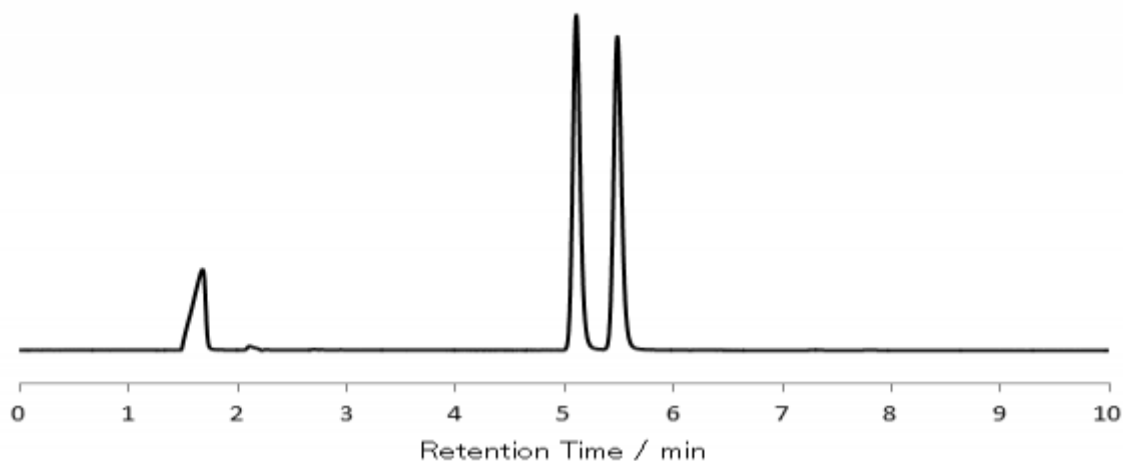

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	water / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.21	Temperature (°C)	25
Injection Volume (µL)	0.0005 (mg)	Sample Concentration	
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.1	1.29	1 - 2	2.91
2	5.5	1.46		