


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Perisoxal

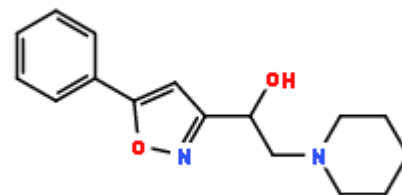

## CHROMATOGRAM

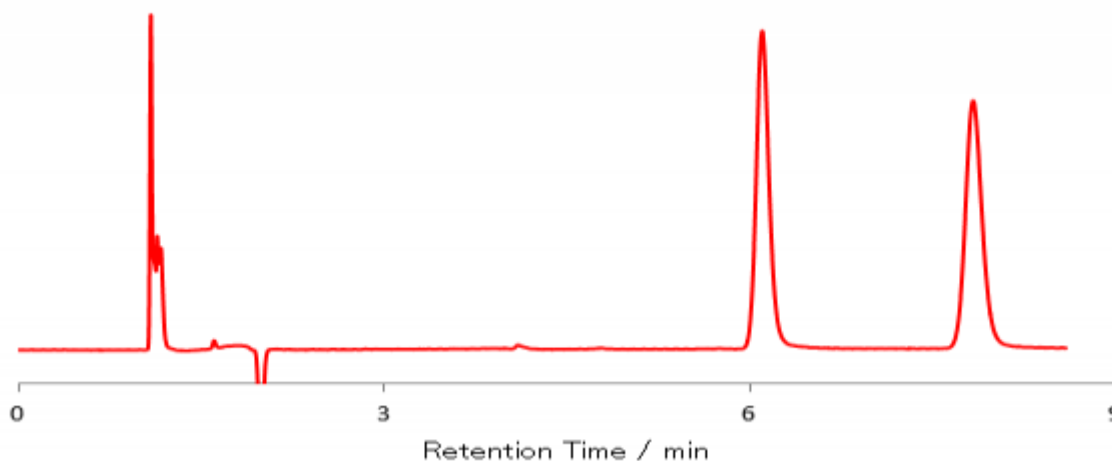

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IC-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	n-hexane / 2-propanol / diethylamine = 90 / 10 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	0.425	Temperature (°C)	25
Injection Volume (µL)	0.0001 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.1	4.69	1 - 2	7.93
2	7.8	6.30		