


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

N-benzoyl-DL-Leucine

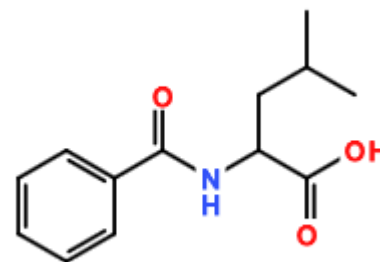

## CHROMATOGRAM

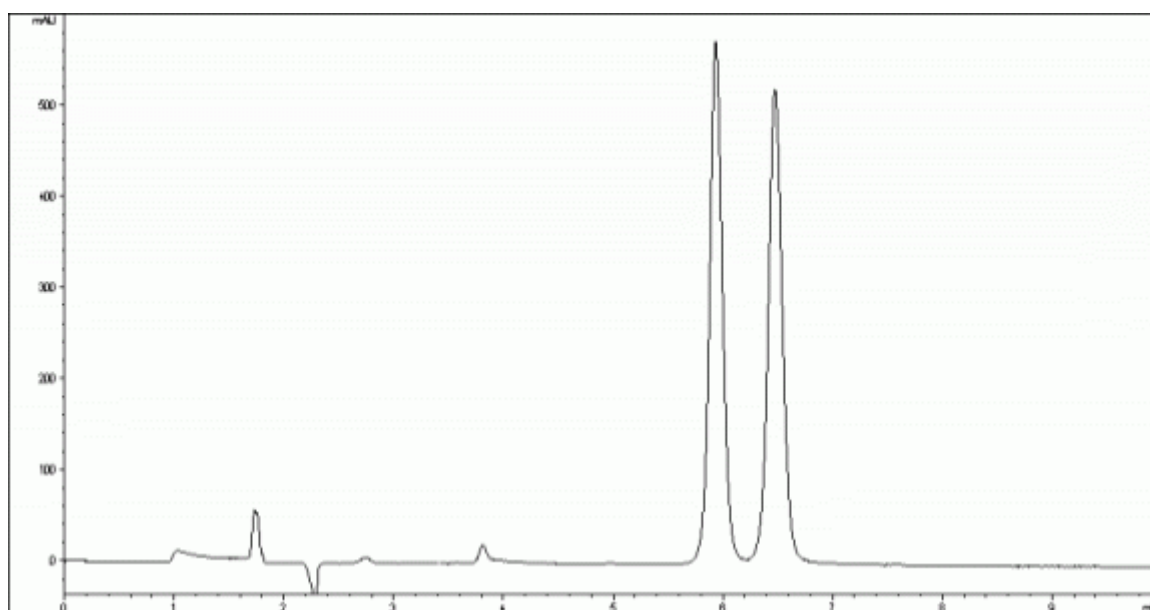

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	n-hexane / dichloromethane / ethanol / trifluoroacetic acid = 50 / 50 / 2 / 0.1		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1	Temperature (°C)	25
Injection Volume (µL)	0.005 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	5.9	2.49	1 - 2	2.47
2	6.5	2.81		