


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

N-CBZ-DL-Valine

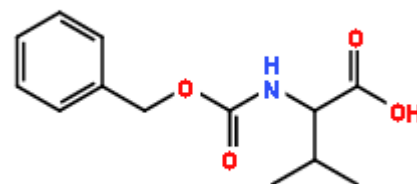

CHROMATOGRAM

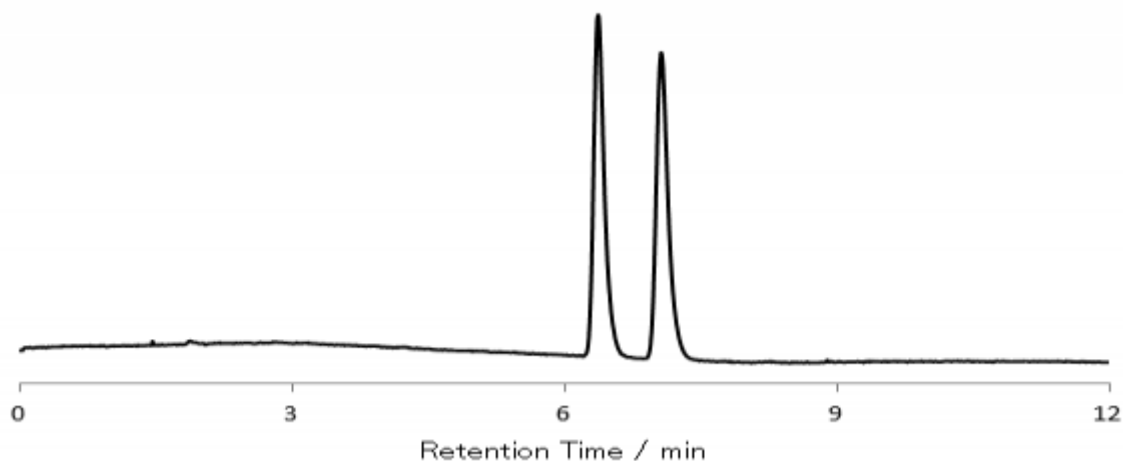

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IH-U		
Column Size (Particle Size)	3 x 100 mm (1.6 µm)		
Mobile Phase Composition	formic acid aq. pH2.0 / acetonitrile = 60 / 40		
Chromatographic Mode	Reverse Phase		
Flow rate (mL/min)	0.425	Temperature (°C)	25
Injection Volume (µL)		Sample Concentration	
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	6.4	1.98	1 - 2	3.15
2	7.1	2.31		