


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Flavanone

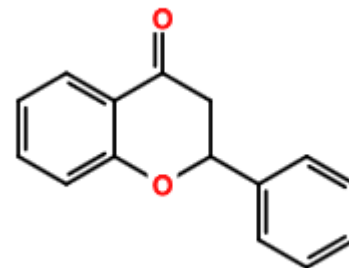

## CHROMATOGRAM

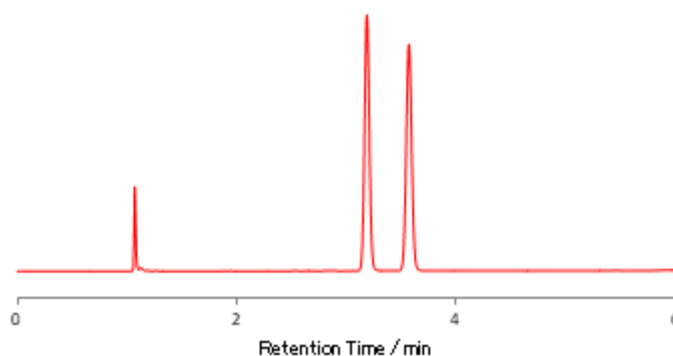

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK <sup>®</sup> IC-U		
Column Size (Particle Size)	3 x 100 mm (1.6 $\mu$ m)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	0.425	Temperature ( $^{\circ}$ C)	25
Injection Volume ( $\mu$ L)	0.0001 (mg)	Sample Concentration	
Detection	UV-VIS 230 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.2	1.97	1 - 2	4.44
2	3.6	2.33		