


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Panthenol

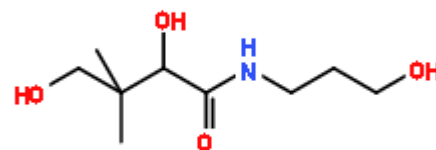

CHROMATOGRAM

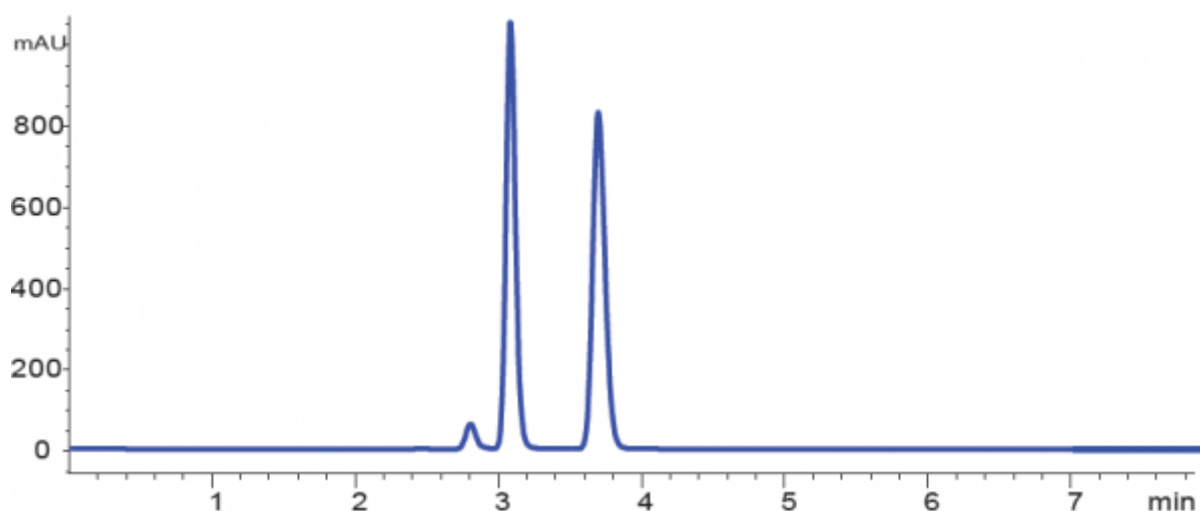

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IG-3		
Column Size (Particle Size)	4.6 x 150 mm (3 µm)		
Mobile Phase Composition	100mM ammonium bicarbonate aq. / methanol = 5 / 95		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.7	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	2.0 mg/mL in Mobile phase
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	3.0		1 - 2	4.30
2	3.6			