


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Flurbiprofen

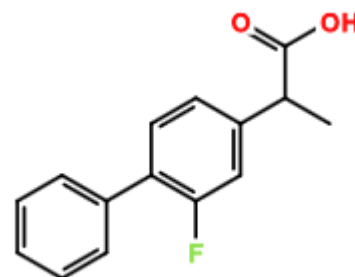

CHROMATOGRAM

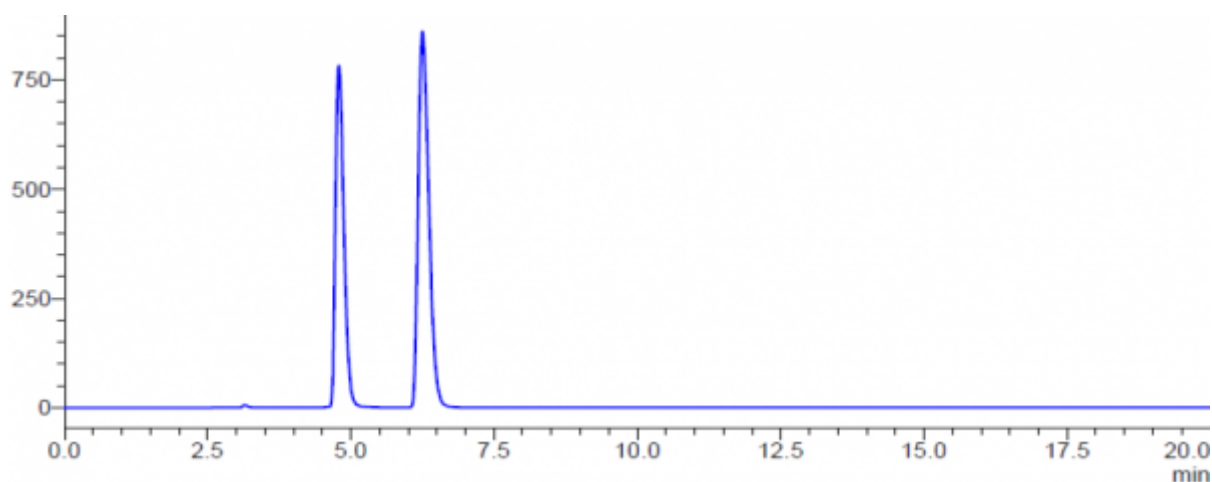

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK [®] IG		
Column Size (Particle Size)	4.6 x 250 mm (5 μm)		
Mobile Phase Composition	n-hexane / ethanol / methanol = 80 / 10 / 10		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (μL)	10	Sample Concentration	0.5 mg/mL in Mobile phase
Detection	UV-VIS 220 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	4.8		1 - 2	4.60
2	6.3			