


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Remdesivir Pentafluoro Phosphoramidite Key Starting Material

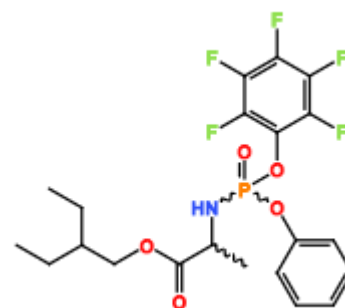

CHROMATOGRAM

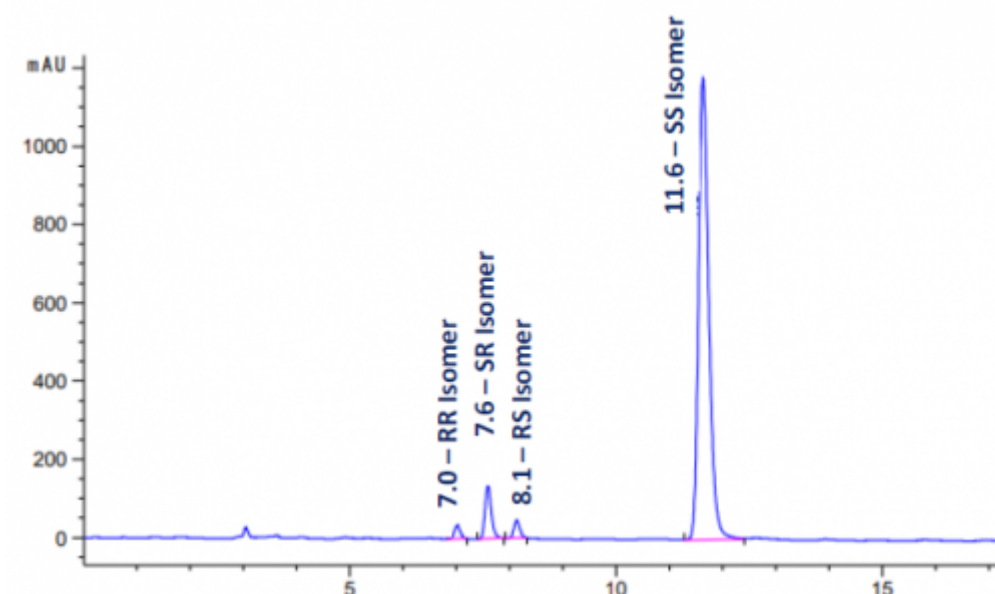

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 90 / 10		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	1 mg/mL in Mobile phase
Detection	UV-VIS 210 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	7.0		1 - 2	2.80
2	7.6		2 - 3	2.60
3	8.1		3 - 4	12.80
4	11.6			