


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

Rotigotine

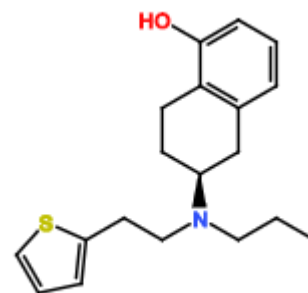

## CHROMATOGRAM

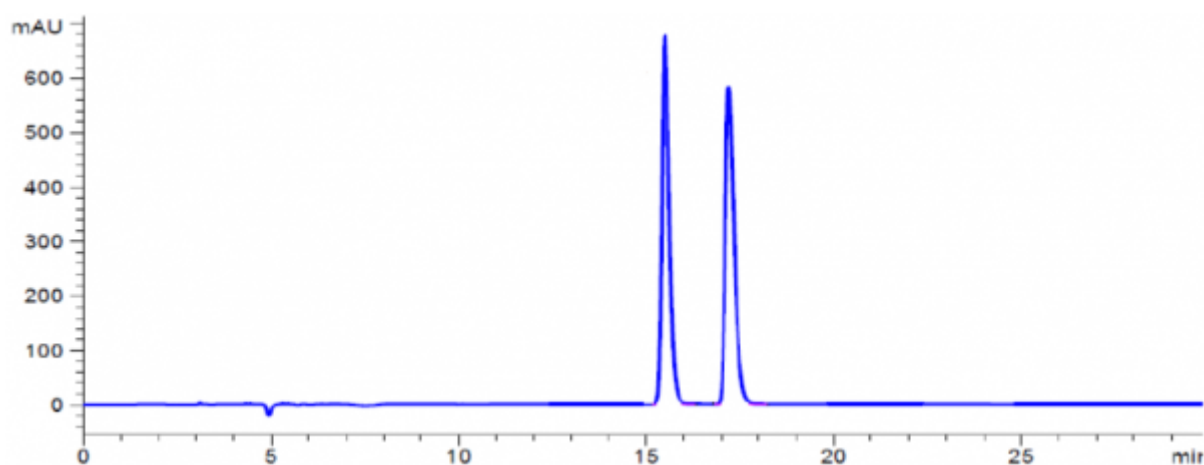

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IB N-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / methanol / diethylamine = 97 / 3 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	10	Sample Concentration	1.0 mg/mL in Mobile phase
Detection	UV-VIS 225 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	15.5		1 - 2	4.30
2	17.2			