


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Remdesivir

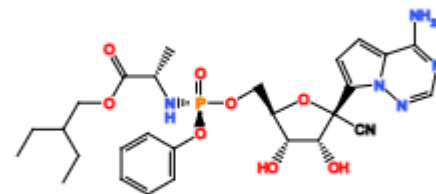

CHROMATOGRAM

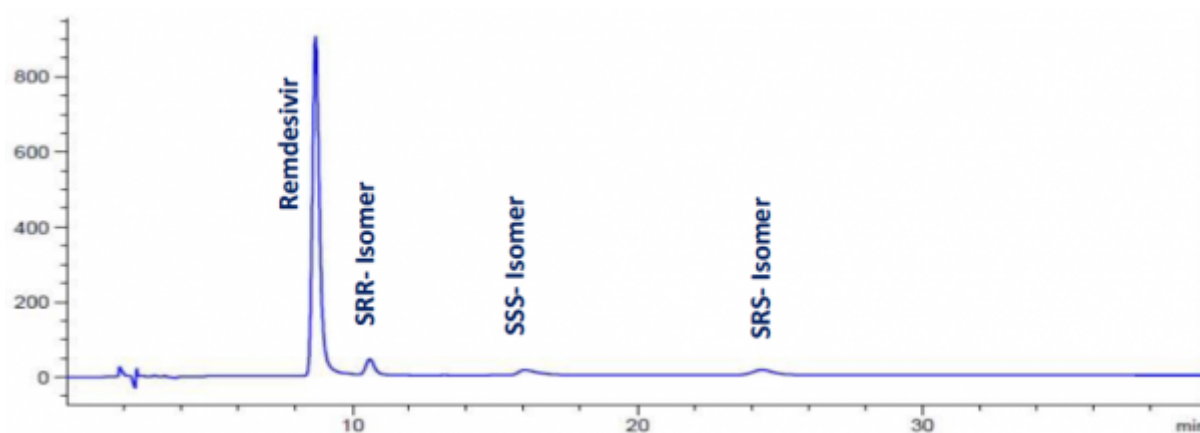

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IA-3		
Column Size (Particle Size)	4.6 x 250 mm (3 µm)		
Mobile Phase Composition	n-hexane / ethanol / 2-propanol / ethanolamine / formic acid = 80 / 5 / 15 / 0.05 / 0.1		
Chromatographic Mode	Normal Phase		

Flow rate (mL/min)	1.5	Temperature (°C)	40
Injection Volume (µL)	10	Sample Concentration	1 mg/mL
Detection	UV-VIS 245 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	8.7		1 - 2	3.90
2	10.6		2 - 3	6.20
3	16.1		3 - 4	6.60
4	24.4			

Comment

Analyte was dissolved in n-hexane/ethanol (1:1).