


APPLICATION RESULTS


COMPOUND IDENTIFICATION

Compound Name

Netarsudil

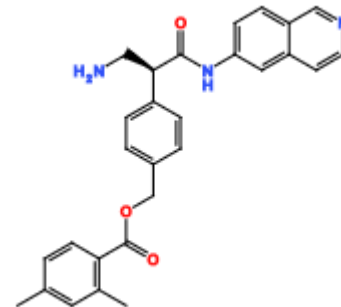

CHROMATOGRAM

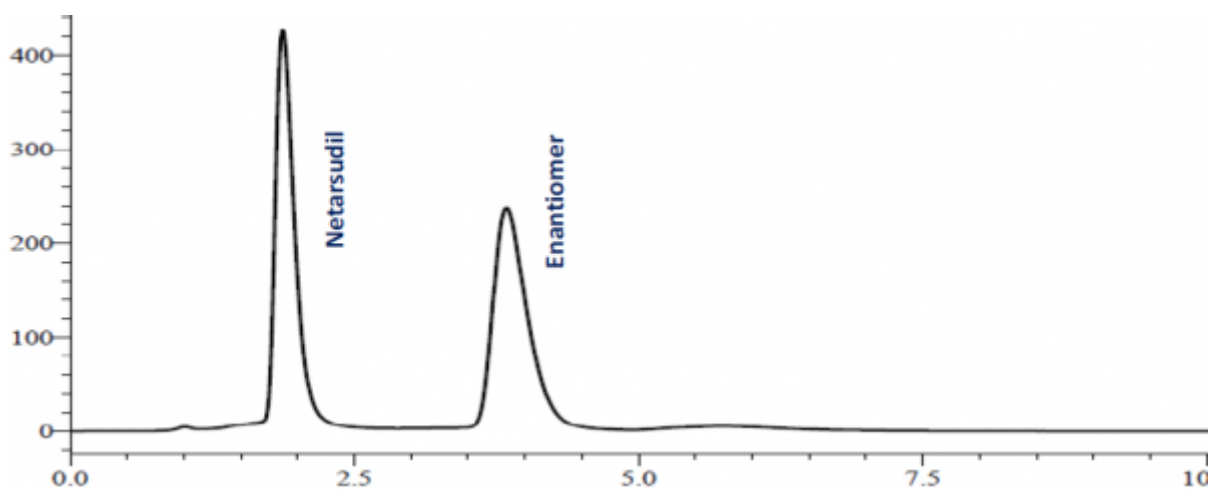

ANALYTICAL REPORT

Column Trade Name	CHIRALPAK IG-U		
Column Size (Particle Size)	3 x 50 mm (1.6 µm)		
Mobile Phase Composition	methanol / acetonitrile / diethylamine = 50 / 50 / 0.1		
Chromatographic Mode	Polar Phase		

Flow rate (mL/min)	0.3	Temperature (°C)	25
Injection Volume (µL)	3	Sample Concentration	0.1 mg/mL in Mobile phase
Detection	UV-VIS 254 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	1.9		1 - 2	4.60
2	3.8			