


## APPLICATION RESULTS


## COMPOUND IDENTIFICATION

Compound Name

gamma-Phenyl-gamma-butyrolactone

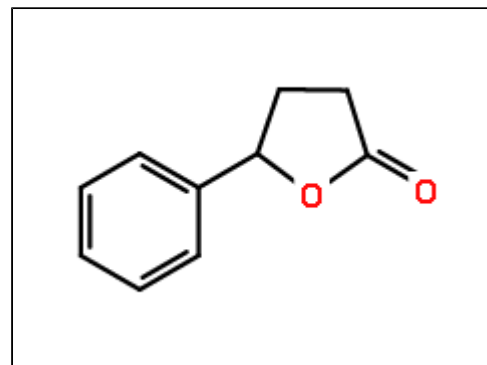

## CHROMATOGRAM

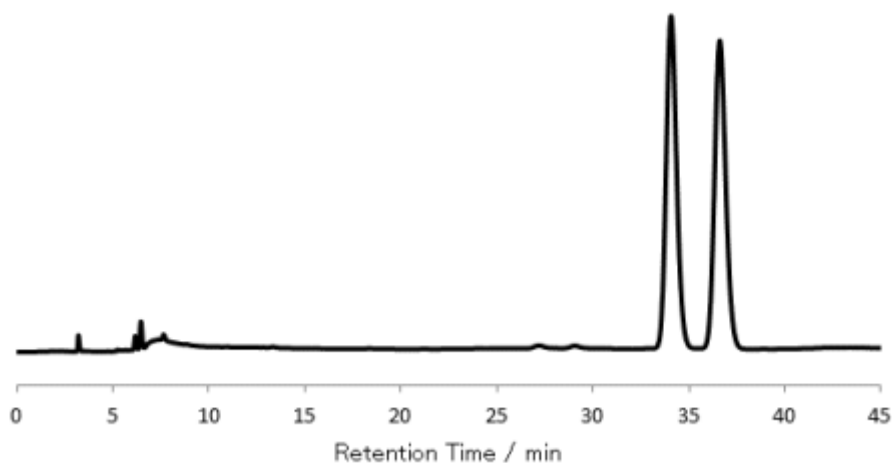

## ANALYTICAL REPORT

Column Trade Name	CHIRALPAK® IK		
Column Size (Particle Size)	4.6 x 250 mm (5 µm)		
Mobile Phase Composition	n-hexane / 2-propanol = 93 / 7		
Chromatographic Mode	Normal Phase		
Flow rate (mL/min)	1.0	Temperature (°C)	25
Injection Volume (µL)	2.0	Sample Concentration	1.5 mg/mL in Mobile phase
Detection	UV-VIS 204 nm		
Elution Order			

No.	Rt	K'	Alpha Value	Resolution
1	34.1	10.36	1 - 2	2.43
2	36.6	11.21		